



MILKOTRONIC LTD


www.lactoscan.com
www.milkotronic.com

LACTOSCAN LW/LAW

MILK ANALYZER Directly controlled by MS Windows

Quick Start Guide

1 UNPACKING INSTRUCTIONS

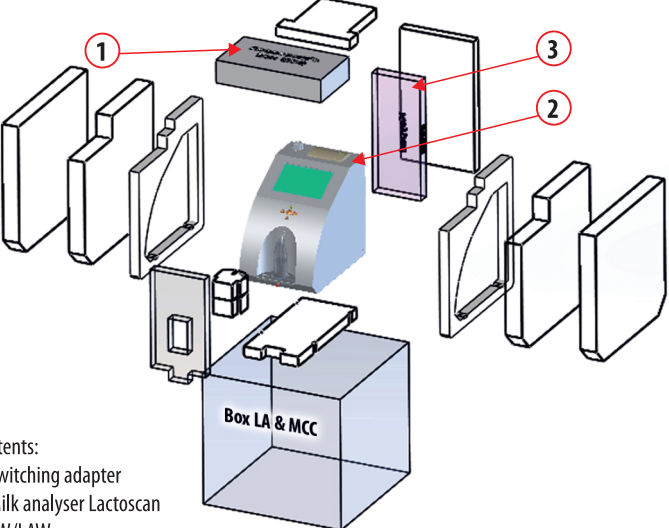

Contents:

1. Switching adapter
2. Milk analyser Lactoscan LW/LAW
3. Box with accessories

2 SAFETY INSTRUCTIONS

1. Read this manual carefully and make sure that you understand all the instructions.
2. For safety purposes the device is equipped with grounded power cable.
If there is no grounded electrical outlet where the device will be used, please, install such before using the device.
3. Place the device on leveled and stable plate.
In case it falls or is severely shocked it may be damaged.
4. Connect to the electrical network in such a way that the power cable to stay away from the side for accessing the device and not to be stepped on.
5. Every time before cleaning the device switch it off and unplug it from the electrical outlet.
The device has to remain unplugged till the cleaning completion.
6. Do not disassemble the unit in order to avoid possible electrical shock.
In case of malfunction contact your local dealer.
7. Handle the liquids the device works with carefully, following all the instructions for their preparation.
8. Place the switching adaptor in such a way as to be protected from overflow and spillage of liquids.

3 ASSEMBLING THE DEVICE

Open the box. Carefully take out the contents of the box.
Place the machine on flat level surface





4 SWITCHING ON THE SYSTEM

1. Switch on the usb hub on the back of the device (separate power supply)(if available)




Automatic cleaning (option)

2. Plug in the analyser power supply in one of the back panel inlets
3. Plug in any USB device which will be needed (keyboard, mouse, etc.) in the USB connector
4. Switch on the Lactoscan from the back panel switch
5. Power On the Tablet (the left button on the long top side of the tablet)
6. Wait until the software opens


7. If everything is connected properly you should receive a "Ready to Start" message on the screen

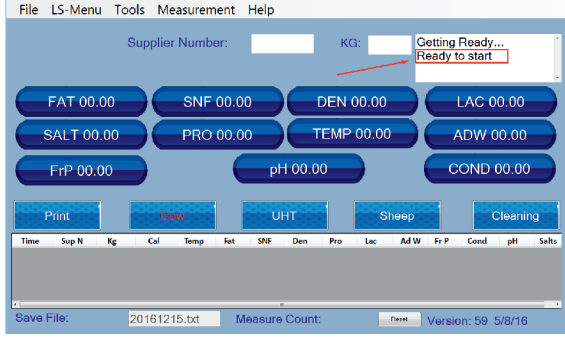


8. Fill the glass with a milk sample and place it in the analyser, when the long hose is in the milk, the short hose is not covered in milk and the glass is placed on the white holder



9. Start a measurement by pressing the calibration name on the screen for the type of milk you are going to measure (if you are measuring cow milk, press the Cow Button)

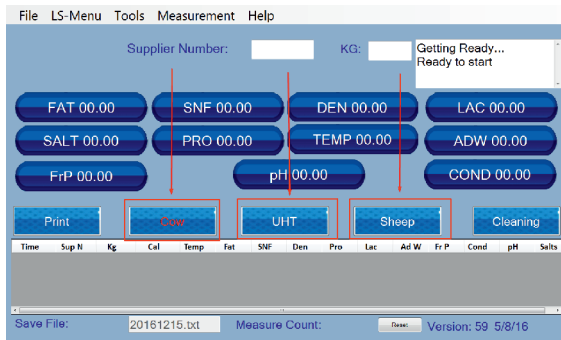


10. When the measurement is started on the screen appears

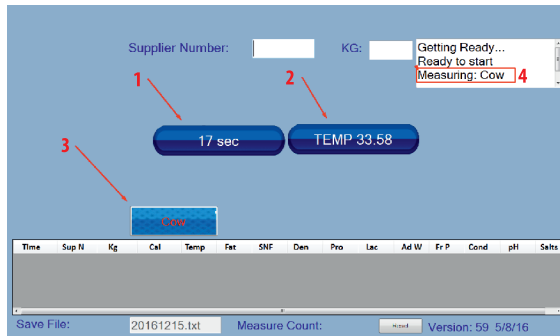

- Where:
- Position 1: The seconds left until the measurement is finished
 - Position 2: Temperature of the sample
 - Position 3: The calibration channel that is being used
 - Position 4: Message about what calibration you are using for measurement


11. When the measurement is finished

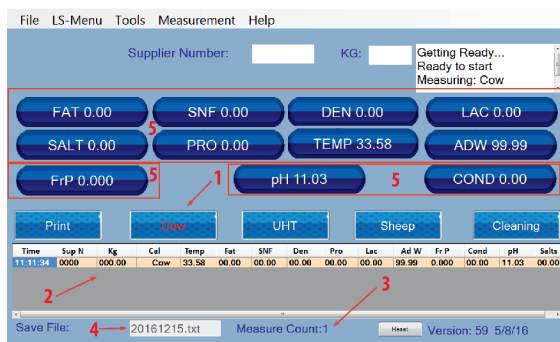

- Where:
- Position 1: Highlighted the last used calibration
 - Position 2: Result collector for all the measurements done since the start of the program
 - Position 3: Measure count since the start of the program
 - Position 4: Archive file for the results
 - Position 5: Fields for the parameter results of the last measurement

12. Continue work by placing a new sample in the glass and press the calibration buttons

13. When finished work put clean water in the glass and start cleaning by pressing the cleaning button

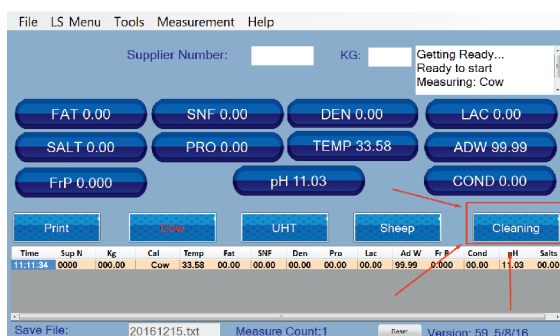


14. Please refer to the user manual about the proper way of cleaning the machine

IMPORTANT


THE MAIN REASON FOR MALFUNCTIONING OF THE DEVICE IS THE BAD CLEANING OF THE SYSTEM AFTER MAKING ANALYSIS. In case of malfunction due to the bad cleaning of the analyser your guarantee is not valid anymore and any repair has to be paid.

5 TURNING OFF THE SYSTEM

1. To turn off the tablet select **File -> Shutdown**

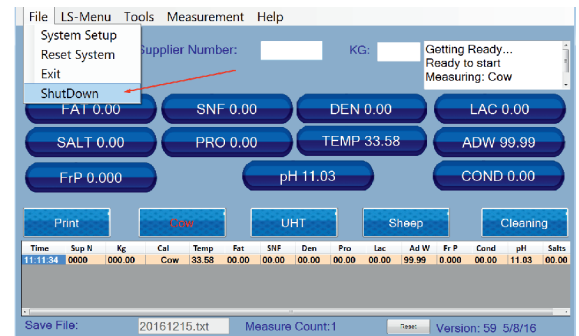

2. To exit the program select **File->Exit** to close the software to use the tablet

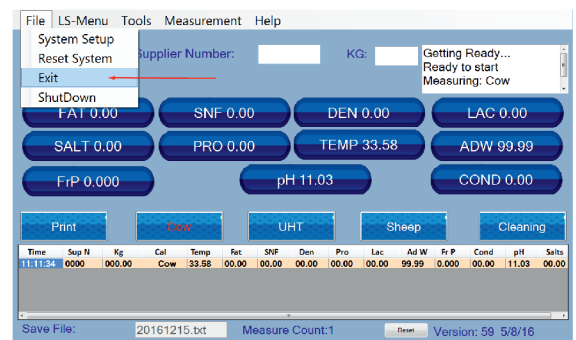

3. Turn off the analyser from the back panel switch


*NOTE: If you have any questions about any of the procedures, please contact Milkotronic LTD

6 OPTIONS AND ACCESSORIES



7 MORE INFORMATION

You can download the Instruction Manual of the device from our website http://www.lactoscan.com/editor/info/manuals/LW/Lactoscan_LW_EN.pdf
You can also find it in the Help Menu of the analyser:

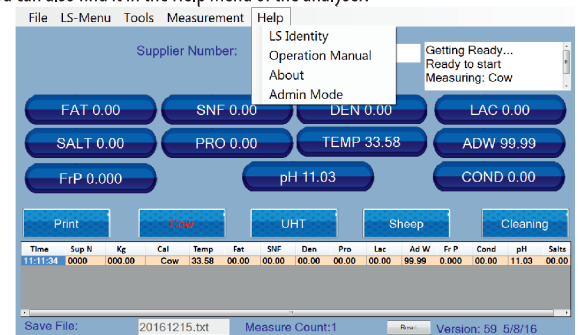

The Instruction Manual of the Data Base Application can be downloaded from the Help Menu of the analyser.

8 SUPPORT

For any additional information and further question regarding the servicing of the device, contact us at: mainservice@lactoscan.com

HEAD OFFICE:
4, Narodni Buditeli Str.
8900 Nova Zagora
BULGARIA
Phone/Fax: + 359 457 67082

MILKOTRONIC LTD
office@lactoscan.com
www.lactoscan.com
www.milkotronic.com