

Milkotronic Ltd

LACTOSCAN LA

MILK ANALYZER

LCD display – 4 lines x 16 characters

Operation manual

Switching Adapter

- **Input:** 100 - 240 V ~1.6 A max.
50-60 Hz
- **Output:** +12 V \equiv 4.17A min.
- **Output power:** 50 - 65 W

Measurement modes

- | | |
|----------------------------|-------------------------------------|
| • cow milk | <input checked="" type="checkbox"/> |
| • sheep milk | <input checked="" type="checkbox"/> |
| • UHT milk | <input checked="" type="checkbox"/> |
| • goat milk | <input type="checkbox"/> |
| • buffalo milk | <input type="checkbox"/> |
| • cream | <input type="checkbox"/> |
| • whey | <input type="checkbox"/> |
| • recovered milk | <input type="checkbox"/> |
| • other /pasteurized milk/ | <input type="checkbox"/> |

CAUTION!

Keep the switching adapter dry!
Please, read and follow strictly all the instructions in the manual.

Due to continuous improvement in milk analyser, information contained in this manual is subject to change without notice. Contact producer for revisions and corrections

4, Narodni Buditeli Str.
8900 Nova Zagora
BULGARIA
Phone/Fax: + 359 457 67082
e-mail: office@lactoscan.com
www.lactoscan.com
www.milkotronic.com

SAFETY INSTRUCTIONS

- 1. Read this manual carefully and make sure that you understand all the instructions.**
- 2. For safety purposes the device is equipped with grounded power cable. If there is no grounded electrical outlet where the device will be used, please, install such before using the device.**
- 3. Place the device on leveled and stable plate. In case it falls or is severely shocked it may be damaged.**
- 4. Connect to the electrical network in such a way that the power cable to stay away from the side for accessing the device and not to be stepped on.**
- 5. Every time before cleaning the device switch it off and unplug it from the electrical outlet. The device has to remain unplugged till the cleaning completion.**
- 6. Do not disassemble the unit in order to avoid possible electrical shock. In case of malfunction contact your local dealer.**
- 7. Handle the liquids the device works with carefully, following all the instructions for their preparation.**
- 8. Place the switching adaptor in such a way as to be protected from overflow and spillage of liquids.**



PARTS AND ACCESSORIES

In the table below the standard delivery configuration of the milk analyser is listed:

No	Description	Item №	pcs
1.	Ultrasonic portable milk analyzer	LSLA001	1
1 sample measurement time		50 sec	<input type="checkbox"/>
		30 sec	<input type="checkbox"/>
2.	Operation manual		
3.	Plastic sample holder	LSLA003	1
4.	Spare Pipes	LSLA004	2
5.	12 V DC Power Supply Cable	LSLA005	2
6.	Alkaline cleaning solution Lactodaily	100 g	1
7.	Acidic cleaning solution Lactoweekly	100 g	1

In the table below the milk analyser spares and accessories, which are delivered on customers' request are listed:

No	Description a) included in the set: <input checked="" type="checkbox"/> b) not included in the set (may be additionally bought): <input type="checkbox"/>	Item №	pcs	<input checked="" type="checkbox"/> / <input type="checkbox"/>
8.	RS232 Interface Cable - Analyser-IBM PC	LSLA006		<input type="checkbox"/>
9.	Service Pack - CD	LSLA007		<input type="checkbox"/>
10.	pH measuring system	LSLAS009	1	<input type="checkbox"/>
11.	pH probe with cable and holder	LSLA010	1	<input type="checkbox"/>
12.	Buffer solution Ph 60 ml (pH7.00±0.01/20°C)	LSLA011	1	<input type="checkbox"/>
13.	Buffer solution pH 60 ml (pH4.00±0.01/20°C)	LSLA012	1	<input type="checkbox"/>
14.	Milk conductivity measuring system	LSLA013	1	<input type="checkbox"/>
15.	Buffer solution conductivity 50 ml (5.02 (±5%) mS/cm (18±0.1°C)	LSLA014	1	<input type="checkbox"/>
16.	Real time clock	LSLA015	1	<input type="checkbox"/>
17.	High-fat measuring function	LSLA016	1	<input type="checkbox"/>

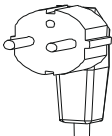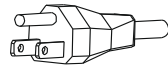18.	ECS POS Serial Printer	LSLA017	1	<input type="checkbox"/>
19.	12 V Serial Printer Power Supply Cable	LSLA018	1	<input type="checkbox"/>
20.	RS232 Interface Cable - Milk Analyser – Serial Printer	LSLA019	1	<input type="checkbox"/>
21.	Power supply from embedded accumulator with charger	LSLA021	1	<input type="checkbox"/>
22.	Plug type		1	<input checked="" type="checkbox"/>
			1	<input type="checkbox"/>
23.	Spare O-ring for the pH probe		1	<input type="checkbox"/>
24.	USB flash drive option		1	<input type="checkbox"/>

1. FUNCTION

The function of the milk analyser is to make quick analysis of milk on fat (FAT), non-fat solids (SNF), proteins, lactose and water content percentages, temperature (°C), pH, freezing point, salts, total solids, conductivity as well as density of one and the same sample directly after milking, at collecting and during processing.

2. TECHNICAL PARAMETERS

2.1. Working modes characteristics:

The program of the milk analyser has four working modes.

2.1.1. Measurement mode milk / dairy product – first type

2.1.2. Measurement mode milk / dairy product – second type

2.1.3. Measurement mode milk / dairy product – third type

These modes have been calibrated on customers' request for 3 milk types from the following: cow, sheep, UHT, buffalo, goat, camel milk, cream, ice cream mixtures, whey, recovered milk, etc. before leaving the production facilities and the text on the display will be for the corresponding types, as is indicated on page 2 Measurement modes.

2.1.4. Cleaning

2.1.4.1. Current

2.1.4.2. Final

Note:

For devices with a 5 button keypad on the front panel by pressing the extra buttons, you can run the following commands: Pressing the button labeled A - start re-printing the results Pressing the button labeled B - start of flushing (2.1.4.) These commands are possible when the unit is in Idle mode, i.e. expecting command for measuring sample.

2.2. Measuring range:

Fatfrom 0.01% to 25% (45%*)
SNFfrom 3% to 40%
Density **from 1000 to 1160 kg/m ³
Proteinfrom 2% to 15%
Lactosefrom 0.01 % to 20 %
Water contentfrom 0 % to 70 %
Temperature of milk	from 5oC to 40oC (if measurement is 30 sec, then t° is from 15 to 40°)
Freezing point ****from – 0,400 to – 0,700oC
Saltsfrom 0,4 to 4%
PH***from 0 to 14
Conductivity***from 3 to 14 [mS/cm]
Total Solids *from 0 to 50 %

* Option, on customers' request

** Density data are shown in an abbreviated form. For example 27.3 have to be understood as 1027.3 kg/m³. To determine the milk density, write down the result from the display and add 1000.

Example: result 21,20; density = 1000 + 21,20 = 1021,2 kg/m³

The abbreviated form of the density is used also when entering data for samples in working mode **Recalibrate**, for example:

If the measured sample density is 1034.5 kg/m³, then in the menu for entering the samples parameters used for calibration, across the parameter Den = , you have to enter 34.5.

*** pH and conductivity measurements are optional and are embedded in the device on customers' request.

**** Please, carefully read Appendix Freezing Point.

2.3. Accuracy:

Fat $\pm 0.06\%$
SNF $\pm 0.15\%$
Density $\pm 0.3 \text{ kg/m}^3$
Proteins $\pm 0.15\%$
Lactose $\pm 0.20\%$
Water content $\pm 3.0\%$
Temperature of milk $\pm 1^\circ\text{C}$
Freezing point..... $\pm 0.005^\circ\text{C}$
Salts $\pm 0.05\%$
PH ± 0.05
Conductivity ± 0.05
Total solids $\pm 0.17\%$



Accuracy values in point 2.3 are in dependence on the correctness of the corresponding chemical method, used for component content determination. In point 2.3. are used the following reference methods: Gerber – for fat, gravimetric – for SNF, Kjeldahl – for protein. The boundary for maximum variation of repeatability when the power supply voltage is from +10 to – 15% from the nominal voltage values (220V) have to be no more than 0.8 accuracy according point 2.3. The analyser is used in conditions free of outer electrical and magnetic fields (except the magnetic field of the Earth) and vibrations.

2.4 Correct ambient conditions:

Accuracy is guaranteed in case of normal ambient conditions:

Air temperature.....from 10°C to 40°C (43°C)
Relative humidityfrom 30% to 80%
Power supply220V (110V)
extent of contamination at normal environmental conditions.....2

2.5. Dimensions:

.....250/290/300 mm, mass 4,0 kg


2.6. Continuous working time:

.....non-stop

2.7 Milk sample volume per one measurement:

..... 15 cm^3 (= 25 ml)


Fig.1 Front panel

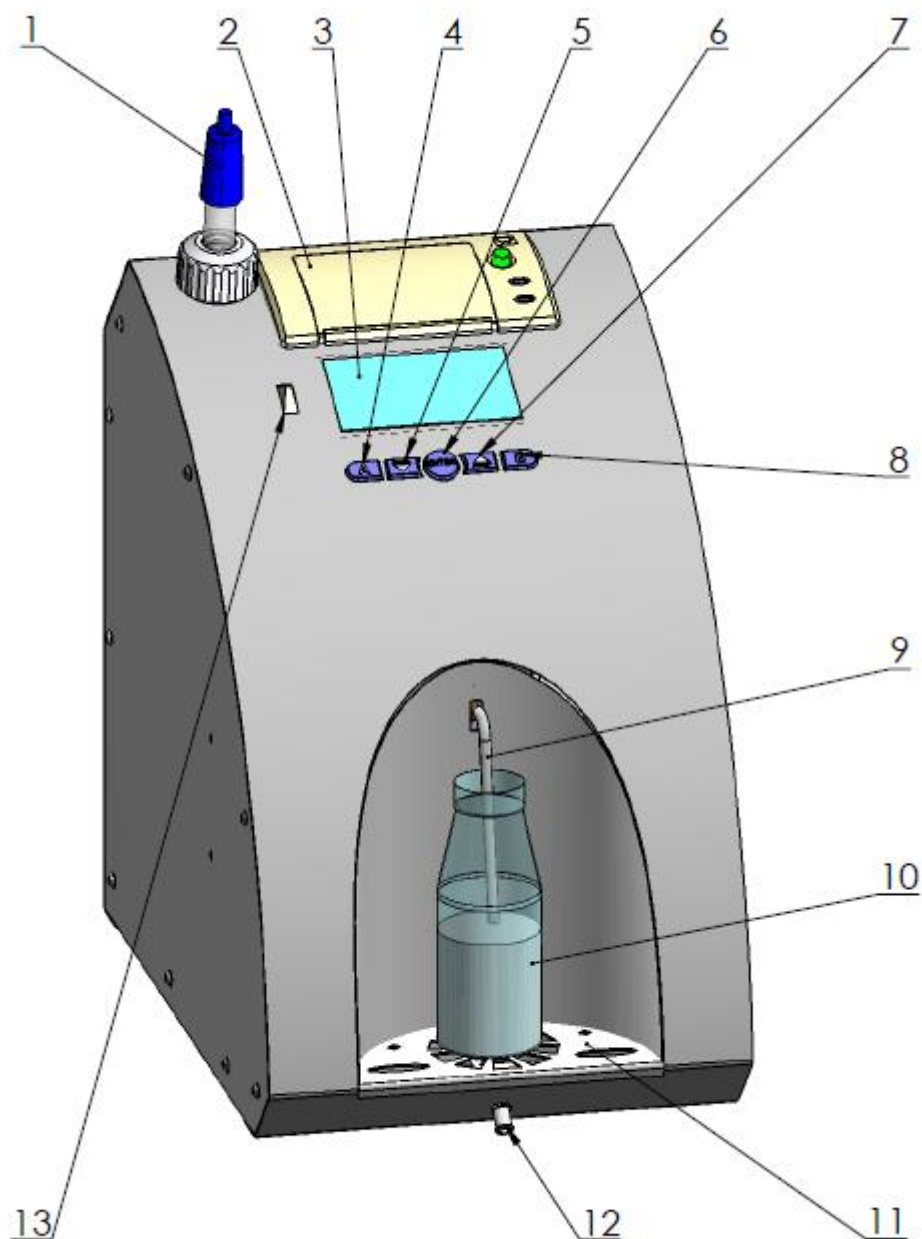

- 1. pH probe & holder (option)
- 2. Printer
- 3. Display
- 4-8 Buttons
- 9. input pipe

- 10. sample holder
- 11. ss grid
- 12. waste liquid pipe
- 13. USB flash-drive slot

Fig.2 Back panel

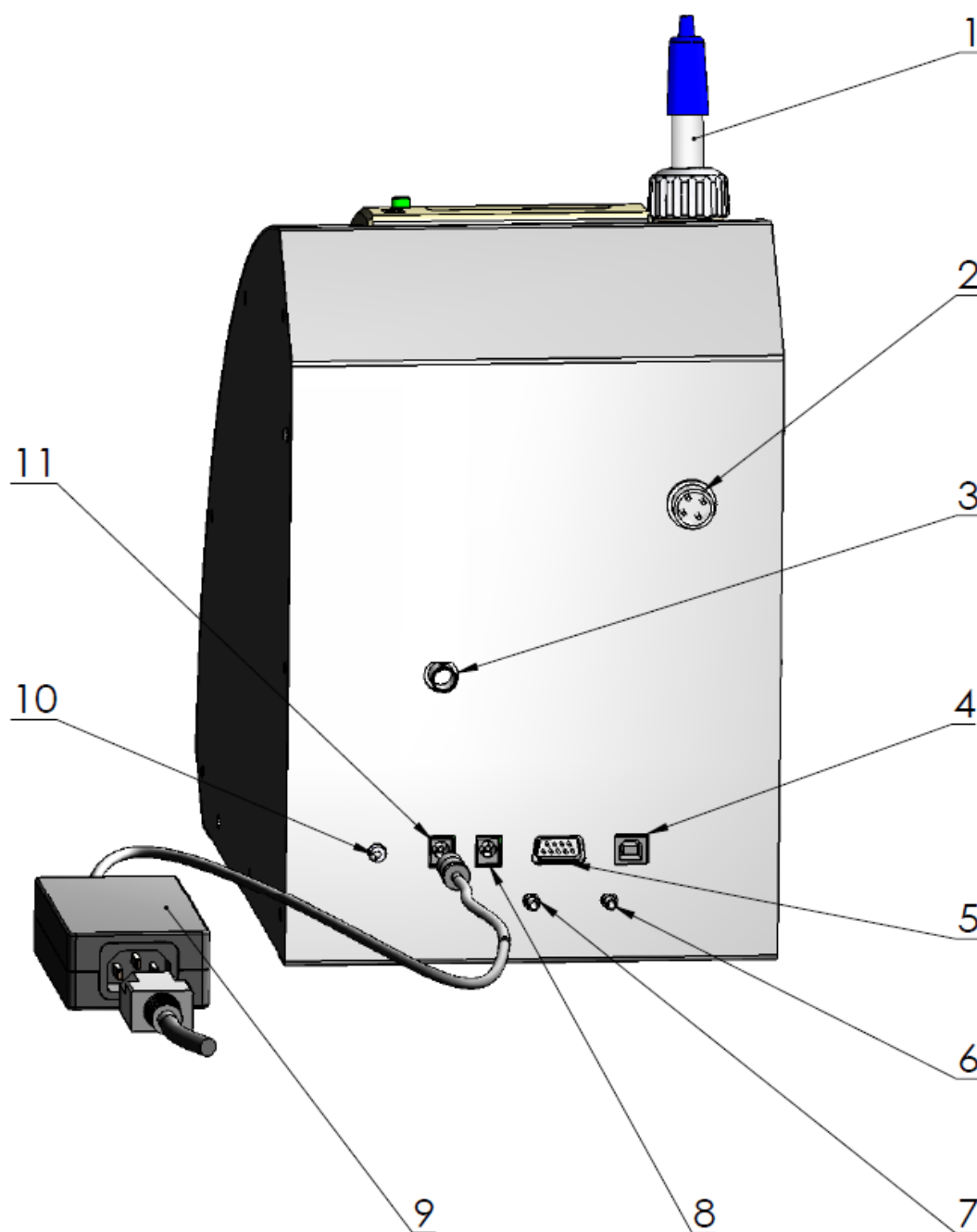

1. pH probe (option)

2. weight scales input (option)

3. pH probe input (option)

4. USB/RS232 (option)

5. RS232/printer

6-7 connectors for the cleaning solution/waste liquid

8.DC Power Supply output

9. Switching adapter

10. Power switch

11. DC Power Supply input

Fig 3 Principle working scheme

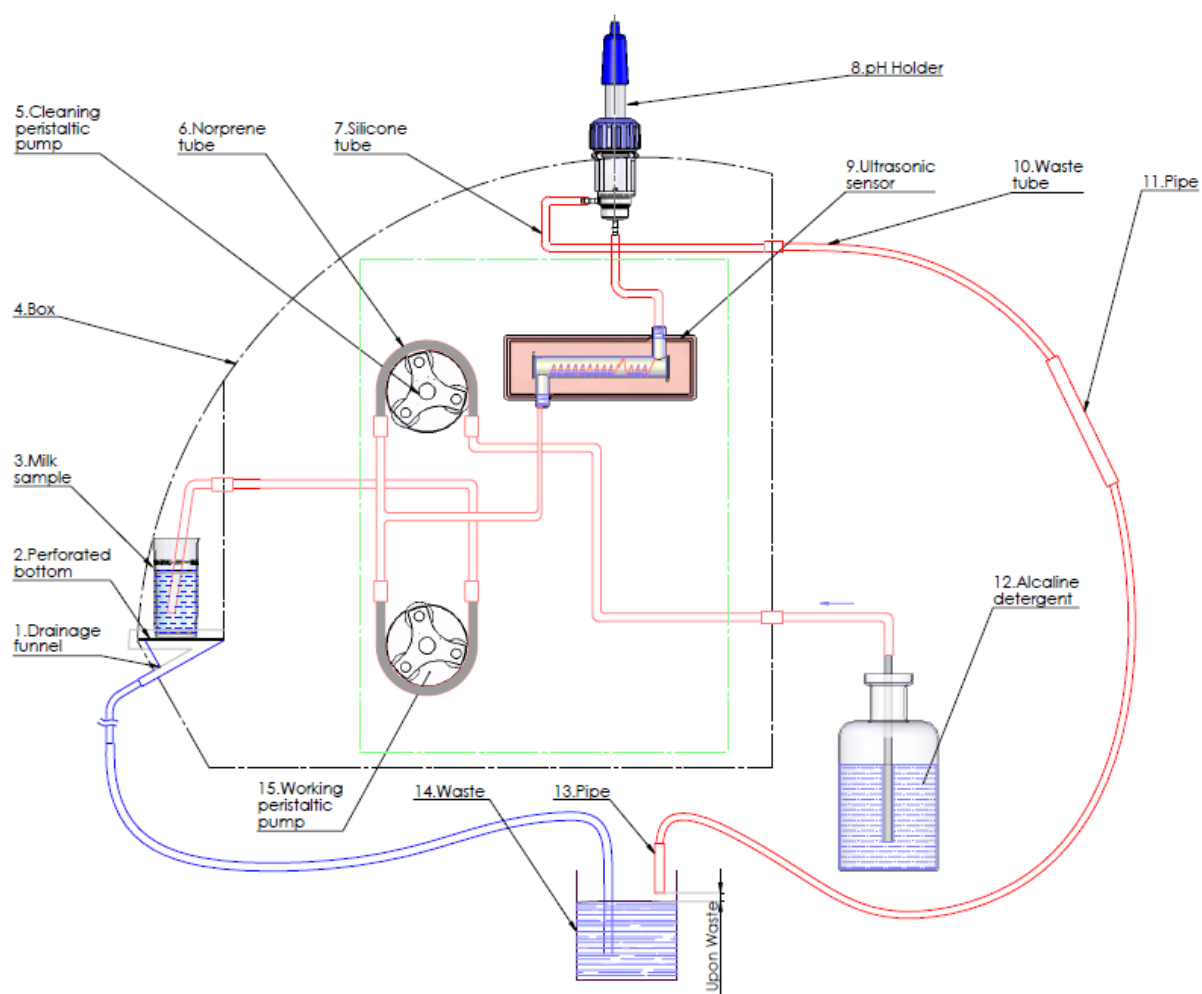


Fig. 4 **Connecting peripheral devices**

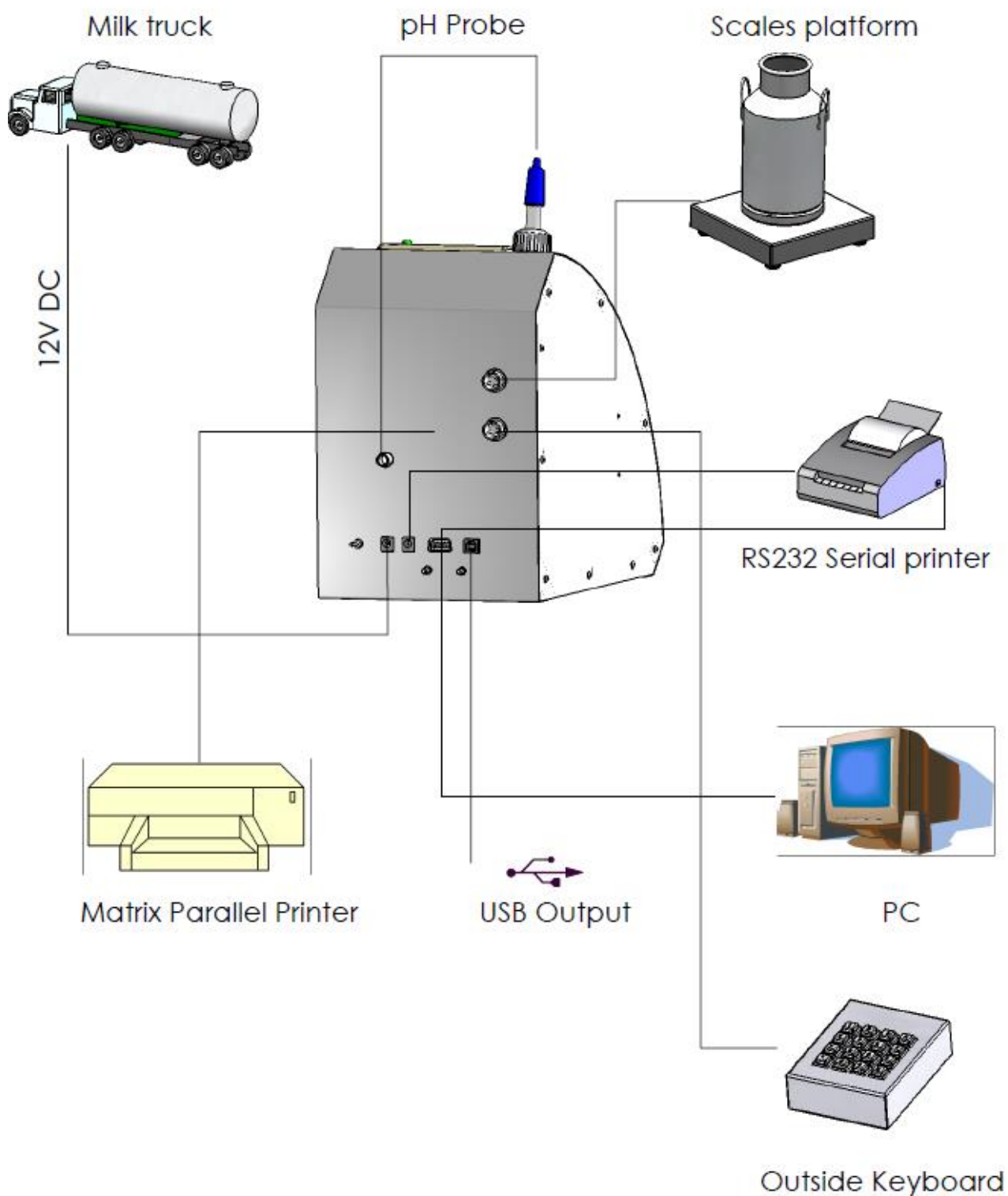
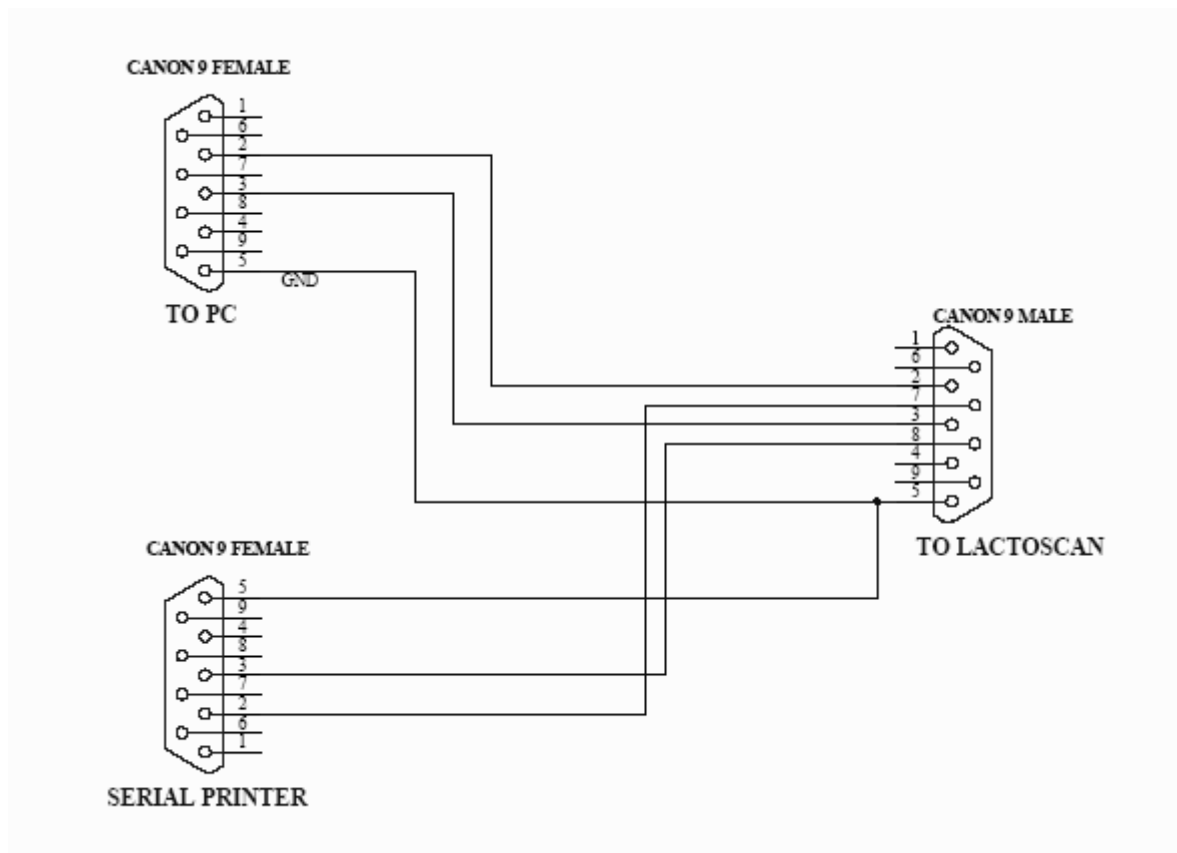


Fig. 5 Cable Description**90-1801-0008**

RS232 Interface Cable - Milk Analyser – Serial Printer/ IBM PC

**90-1801-0009**

DC 12V Power Supply Milk Analyzer Cable

1. GND
2. No connection
3. No connection
4. 12V DC

There's an option in the analyser – embedded USB interface (USB to RS232 Serial Converter) for connection with computer type IBM PC. It is intended for connecting computers from laptop type, which have no RS 232 interface (COM Ports – DB9 Connector). This option may be available together with the standard RS232 – connector DB9, which always exists in the analyzers. At one and the same time, connection analyzer – computer, can be established only through one of the couplings – either only DB9, or only USB, i.e. analyzer could not be connected at one and the same time, using one of the interfaces with one and using the other interface to another, second computer. The coupling for this option (Printer Type) is at the back of the analyzer, next to the standard RS232 DB9 connector. Other USB devices like printers, keypads etc COULD NOT be connected to this coupling.

USB interface is based on the element MCP2200 of the company Microchip Technology Inc. – site: <http://www.microchip.com>. For establishing a connection through this interface a driver has to be installed for MCP2200 in the correspondent computer, which will be connected with the analyzer. Please, follow the procedure, from the Internet site of the producer: <http://www.microchip.com/wwwproducts/devices.aspx?ddocname=en546923>, by choosing the suitable for your operation system driver.

After installing the driver, choose the COM port, which will be used for the real communication. For Win XP it is done by executing the following commands: Start -> Settings -> Control panel -> System -> Hardware -> Device Manager -> Ports (COM and LPT) – Right Click -> Properties -> Port Settings -> Advanced -> COM Port Number. Choose Number 1 or if another number is chosen, for example 3, then, when working with the software tools, in the field COM Port (upper right corner), you have to set the new number of the port, in this case 3.

3. QUALIFICATION OF RAW MILK, THERMALLY TREATED MILK, OTHER DAIRY PRODUCTS AND DERIVATIVES

3.1. Taking samples and preparation for analysis

In order to receive reliable results in qualification of milk, dairy products and derivatives are needed: precise samples taking; correct samples storing (in need to be preserved); correct preparation before making measurement. The rules and requirements for this are described in details in *Appendix Milk sampling*.

3.2. Making the measurement.

3.2.1. Preparing the analyser for working mode

3.2.1.1. Put the analyser on the working place, providing good ventilation and not in the vicinity of heat providing devices or sources. The temperature in the premises has to be in the boundaries 10-40°C.

3.2.1.2. Check if the power switch is in "0" position and that the outlet voltage complies with the voltage indicated on the switching adapter of the analyser

3.2.1.3. Switch on the “**POWER**” button, which starts the identification procedure. For a short time the display shows the number of the software versions, for example:

Milk analyzer xxx
Vers: yy dd-mm-yy
Ser. N. xxxx

Where:

Milk analyser xxx - name of the analyzer, type and time for measurement.

Vers: yy dd-mm-yy – motherboard software version, date for the version

ser. N. xxxx - serial number .

This information is called informational identity (Identity) of the analyser



If in the process of exploitation there is a need to ask a question the company-producer, you have to send the data, written on the display during the above described initialization procedure.

3.2.1.4. Till the analyser is prepared for work (at about 5 minutes) the following message is written on the display: **“Getting ready”**. Above pointed time is in dependence of the environmental temperature and increases with decreasing the temperature.

3.2.1.5. When the device is ready for work, a beep is heard and the display shows: **“Ready to start”**. The analyser is ready to make analysis in mode 1 (normally Cow)

3.2.1.6. If you want to pass to another mode press the button **Enter** and hold it pressed. The following message appears on the display:

**Release button to
start menu**

Release the button **Enter**. The display shows the possible working modes:

Milk selector

Cal1 – Cow
Cal2 – Sheep
Cal3 – UHT

Cleaning
Final clean

Using “up” ▲ and “down” ▼ buttons choose the working mode and press **Enter** in order to start it.

3.2.2. Making analysis



Make sure the analyser is still during measurement. Vibrations may lead to inaccurate results.

To start measurement:

- pour the preliminary prepared sample in the sample holder (now you may use holders with different size - plastic and glass 100 ml);

Fig.6 Input pipe



The input pipe of the analyser is knee-joint (mobile). Move it forward in order to be possible to dip its end in the milk sample and place the sample holder in the recess of the analyser..

- press the button **Enter**.

The analyser sucks the milk, makes the measurement and used milk sample is poured through the holes of the grid below the input pipe of the analyser and enters the vessel for collecting waste liquids (see Fig. 3 Principle working scheme of the analyser). During the measurement the temperature of the sample is shown on the display.

Ignore the results received immediately after switching on the analyser and after measuring distilled water. Make a second measurement with new portion of the same sample.

3.2.3. Displaying the results

3.2.3.1. When the measurement is finished, the sample returns in the sample-holder and the display shows the results. For example:

Results:

F=ff.ff	S=ss.ss
D=dd.dd	P=pp.pp
L=ll.ll	W=ww.ww

Where:

F= ff.ff	- measured FAT in percentage;
S= ss.ss	- measured SNF in percentage;
D= dd.dd	- measured density in percentage;
P= pp.pp	- measured protein in percentage;
L= ll.ll	- measured lactose in percentage;
W= ww.ww	- measured sample's added water in percentage;

By pressing the button "Down" ▼ the display shows the second page, containing the results:

Page 2 Results:
T=tt.tC pH=pp.pp
FP=-0.fff
s=0.sss A=aa.aa

Where:

tt.tC	- sample's temperature;
pp.pp	- sample's pH result – if there is a pH probe connected;
-0.fff	- measured sample's freezing point;
0.sss	- measured salts values;
aa.aa	- measured total solids

By pressing the button "Up" ▲ display shows the third page with results:

Page 3 Results:
L=l.ll

Where:

L= ll.ll	- measured Lactose in %;
-----------------	--------------------------

By pressing the buttons "up" ▲ and "down" ▼, the operator has the possibility to pass from one page result to another.



If the device has an embedded option "Conductivity" and "conductivity measurement" started, the result is shown on the display, showing the basic results replacing lactose results in the following way:

C=xx.xx

In this case the Lactose result is shown on a new page - Page 3 Results.

xx.xx is the measured milk sample's conductivity in [mS/cm]. If the results are outside the limits for this type of sample (see table from the Appendix Conductivity measurement), the cursor flashes after the letter C, reminding that the sample is not correct. On the printout it is printed as !!!.


If the conductivity value is outside measuring range (2-14 mS/cm), the following message appears on the display:

C=OutRg (Out of Range), and on the printout there isn't any line with conductivity value.

3.2.3.2. Write down the results in the form. The results remain on the display till a new measurement is started. If the analyser is connected to a computer, it automatically sends the data to the computer and/or prints them.

4. CLEANING THE ANALYSER

Fig. 7 Connecting the container with detergent

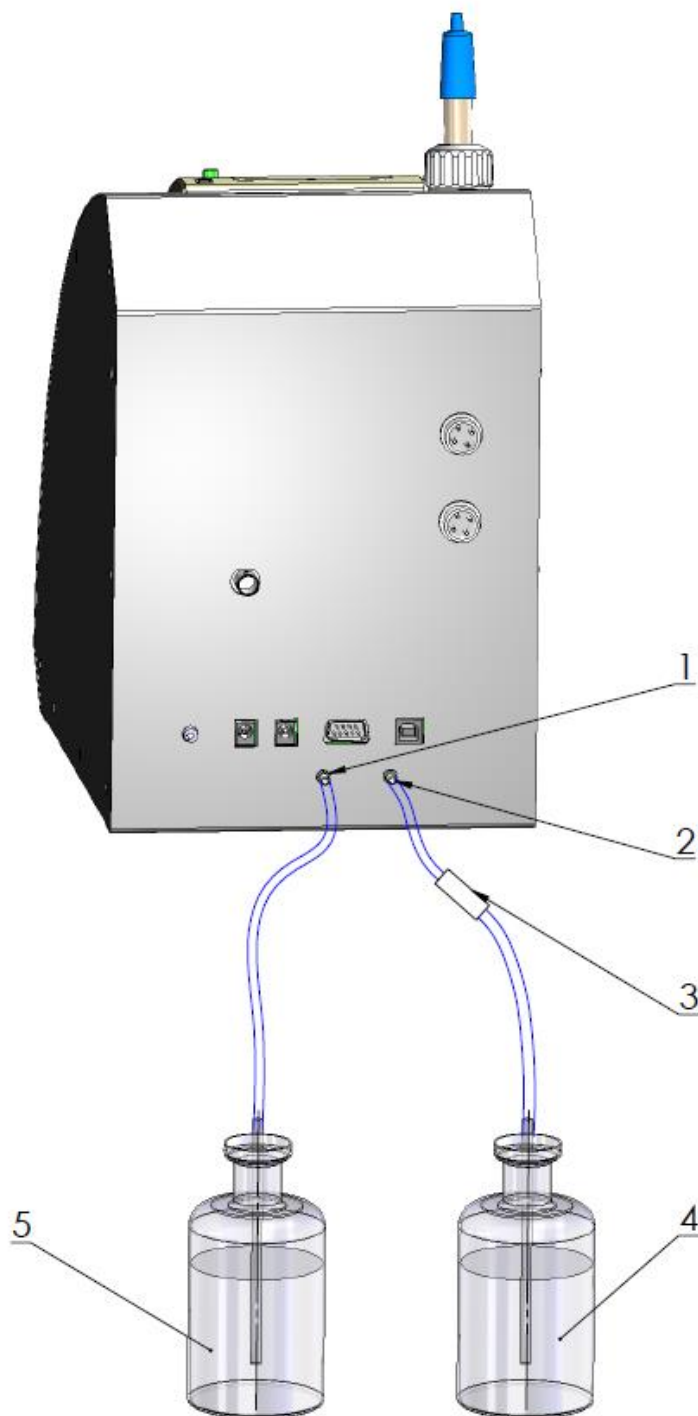

Cleaning solution

Waste solution

- | | |
|----------------------------------------|----------------------------------|
| 1. Connector for the cleaning solution | 4. Vessel for waste solution |
| 2. Connector for the waste liquids | 5. Vessel for the cleaning agent |
| 3. metal pipe | |

This procedure prevents gathering milk fat residues and milk stone on the sensor. The milk stone consists of milk solids, calcium, iron, sulphates, magnesium, etc. All these substances form layer on the pipe and sensor's walls, which leads to deviations in the measurement results and blocking up the piping.



The company-producer recommends usage of the chemicals, supplied with the analyser – alkaline and acidic (Lactodaily and Lactoweekly). They can be ordered separately or together with the analyser. Try to use only these chemicals for cleaning the analyser.

In case you missed to order these chemicals, the alternative is to use alkaline and acidic cleaning solutions for dairy equipment by one the companies, producing such chemicals, as for example:

<http://www.diversey.com>

<http://www.ecolab.com>

<http://www.calvatis.com>

4.1. Automatic cleaning the analyser

Analyser's contamination as a result of the irregular cleaning is the basic reason for inexactness during measurement. In order to be avoid this, in the milk analysers with peristaltic pump is embedded automatic cleaning. For this purpose the analyser has to be prepared on the following way:

The reservoir with the cleaning solution has to be connected towards the analyser and outgoing pipe towards reservoir for collecting the measured samples and used cleaning solution, as is shown on Fig.7.



Pay attention to the pipes in both reservoirs. The pipe in the reservoir with the cleaning solution has to be well dipped in the cleaning solution while the pipe in the reservoir with the already used samples and cleaning solution **MUST NOT** be dipped in the liquid.

Pay attention to the level of the cleaning solution in the reservoir. Add solution when the level is below 2/3!

Description of the cleaning procedures:

There are two ways of cleaning the analyser embedded in it – current and end (final) cleaning.

4.1.1. Current cleaning.

4.1.1.1. Automatic start

The current cleaning is made with alkaline solution (for working solution preparation see 4.2.1.), with which the reservoir for the cleaning solution is filled in.

It is automatically started, without operator's interference after the set time intervals elapse

- 1.** 55 min. after switching on the power supply of the analyser, but idle* work;
- 2.** 15 min. after the last measurement of real milk sample.

*Idle Mode is that part of the standard working mode, when the analyser is not making measurements. There's embedded in the analyser system for measurement of the idle time. The idle time is measured starting from the last action of the operator. In dependence of it (what the operator last did), are taken decisions regarding the cleaning.

There are 2 options:

Option A: If the analyser:

1. Was only switched on but was not started in measurement mode,
2. Or the last action was cleaning,
3. Or the last action was measuring sample with very low Fat (similar to water)

Then the automatic cleaning is started after 55 min.

Option B: If the last thing done with the analyser was measurement of normal milk sample, the automatic cleaning is started after 15 min.

After cleaning a new measurement is started in the above described intervals.

The display shows the following message and a sound signal is emitted:

**Auto clean
started!
Put empty glass**

After this the cleaning is started. The display shows:

**Cleaning
Please wait**

If there is a glass with sample it is completely emptied and then the analyser automatically fulfills procedure for sucking the cleaning solution and rinsing the analyser's inner system. In order the input metal pipe to be cleaned out for a short time a cleaning solution is pumped out for a short time in the already empty glass or in case there is no empty glass there – in the funnel for the waste liquids. It is filled to the middle and then is sucked back and poured in the reservoir with the used liquids. In case that you are sure that you'll not use the analyser for a longer period, but it will be connected to the power supply and switched on, which means that the process of automatic clean will be active, do not leave sample holder in the recess of the analyser, in order not to be overfilled during cleaning from the output pipe. Then the cleaning solution will be directly poured in the funnel for waste liquids.

After the cleaning is finished the displays shows the following:

End of cleaning

After 2 seconds the display shows:

Analyser ready

Then the analyser is ready for normal measurement.

4.1.1.2. Manual start

The current cleaning may be completed by manual start of the menu **Cleaning**. It is used before starting the menu Final Clean. It serves for cleaning the fats from the measuring tract with alkaline cleaning solution with which the canister is filled in. It is started using the menu for choosing the working mode of the analyser. After the measurement is completes, by continuous pressing the button **Enter**, the possible analyser operation modes are shown on the display.

**Cow
Sheep
UHT
Cleaning
Final Clean**

By choosing **Cleaning** the current cleaning is started. The display shows the following:

**Auto clean
started!
Put empty glass**

In this way the current cleaning is started. The procedure and operator's actions are described above.



Do not switch off the device at the end of the working day before the automatic cleaning procedure is completed. If it is not automatically cleaned and there is not a possibility to wait starting the automatic cleaning, then start manually the cleaning procedure with alkaline cleaning solution and if it is necessary with acidic cleaning solution as it is described below.

4.1.2. End (final) cleaning.

4.1.2.1. Final cleaning.

It is done with 3% acidic cleaning detergent (for working solution preparation see 4.2.2.). It is done daily. Serves for cleaning the protein deposits from the measuring system of the analyser, which were not removed during work with the alkaline cleaning solution and removing milk stone layers.



Do not use chemicals, which are not intended for milking equipment and vessels in dairy industry. Pay particular attention to the acidic cleaning solution concentration. The higher concentration may damage the measuring sensor.

Always before the final acidic cleaning rinse the device with alkaline cleaning solution by manual starting if it was not already automatically cleaned. It is necessary to be done in order to remove the milk residues which could react with the acidic cleaning solution.



Do not forget to rinse with clean water!

It is started using the working modes menu. After the measurement is completed, by pressing and holding pressed the button **Enter**, the display shows the possible working modes of the analyser.

**Cow
Sheep
UHT
Cleaning
Final Clean**


1. Choose **Final Clean**. After choosing the final cleaning, the following message appears on the display:

**Put filled with
Water glass
and press Enter
to Continue**

2. Put a glass filled with water.

3. The hose from the Waste pipe has to be placed in the vessel for waste liquids.

Fig 8. Rinsing with water

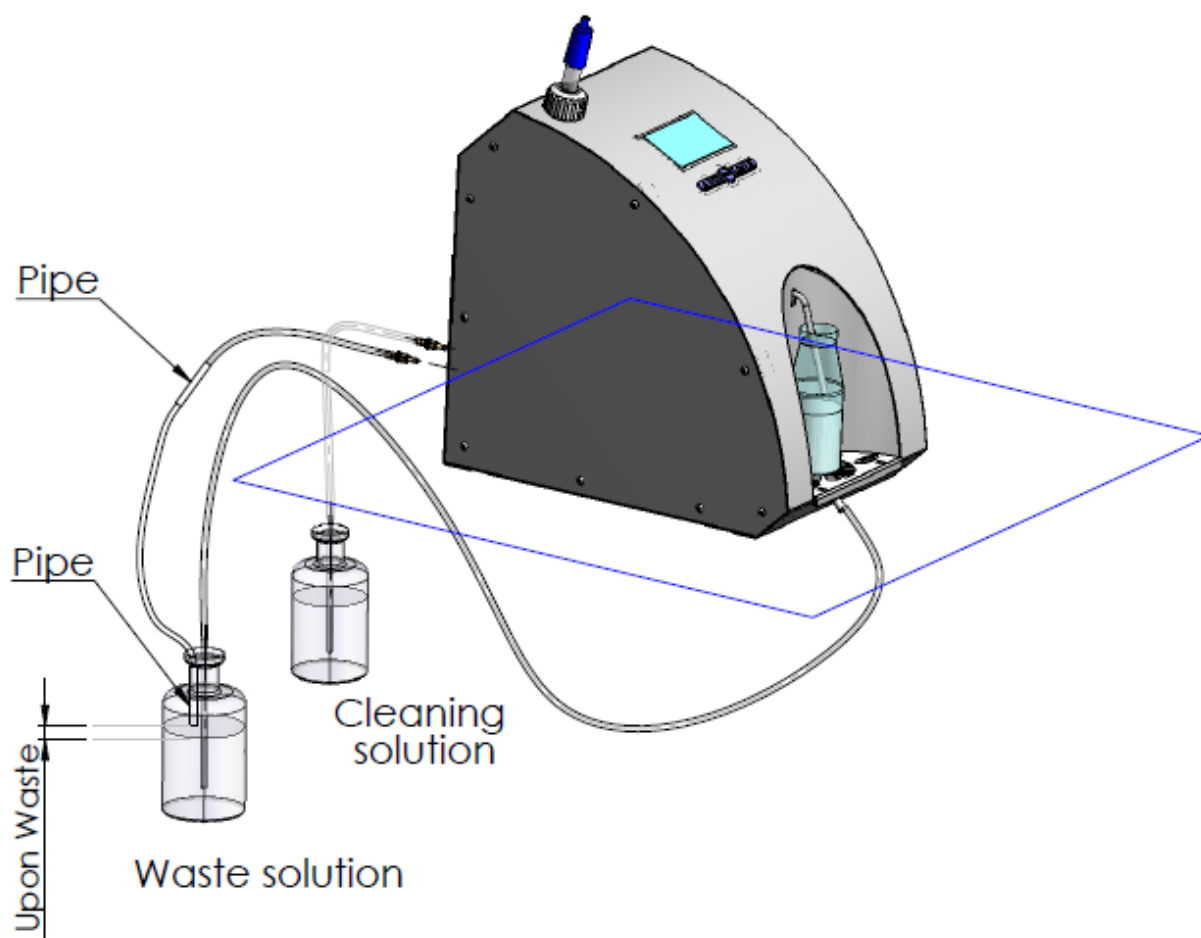

3. Press the button **Enter** in order to continue the procedure. Then the following message appears on the display:

**Cleaning
Please wait**



Do not miss to clean with water because the residues from the alkaline cleaning solution could react with the acidic solution, which will lead to gas and sediment formation.


After rinsing with water, the display shows the following:

**Glass Alcal Deterg
Waste pipe → Glass
press Enter**

5. Place glass filled with alkaline cleaning solution.

6. The pipe, coming from the orifice with inscription Waste Pipe has to be placed in the glass.

Fig 9. Cleaning

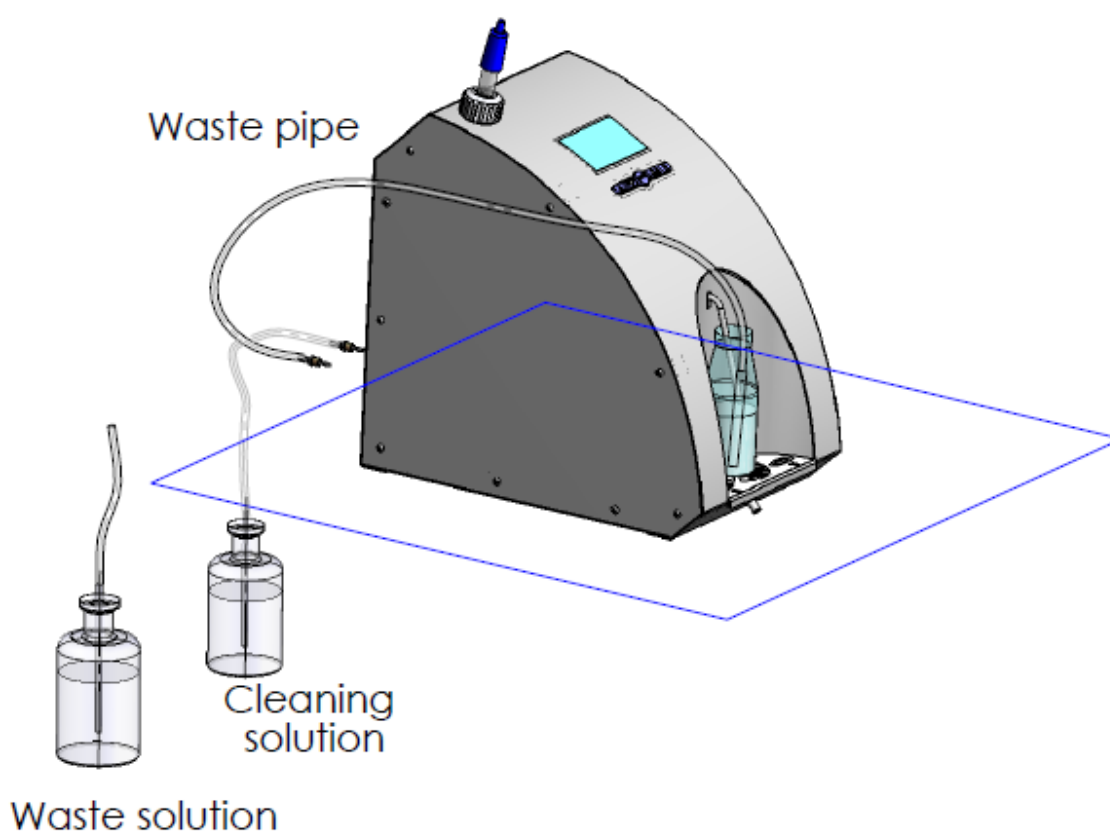

7. Pressing button Enter starts the procedure for cleaning with alkaline solution.

The following appears on the display:

Cleaning: 1

Where, changing number shows number of cycles done by the pump.

8. After the cleaning procedure is finished, the following message appears, reminding the the analyzer needs rinsing.

9. Rinsing procedure is executed, described in p. 2. After it is finished, the display shows:

**Glass Acid Detrg
Waste pipe → Glass
press Enter**

10. Place glass filled with acidic cleaning solution.

11. The pipe, coming from the orifice with inscriptionc Waste Pipe has to be placed in the glass.

12. Pressing button Enter starts the procedure for cleaning with acidic solution.

The following appears on the display:

Cleaning: X

Where, X is changing number showing number of cycles done by the pump.

13. After the cleaning procedure is finished, the following message appears, reminding the analyzer needs rinsing.

14. Rinsing procedure is executed, described in p. 2. After it is finished, the display shows:

**End of Cleaning
Ready to Start**

15. Procedure Final Clean is finished. Now the analyser is ready for normal measurement or to be switched off from the power supply.



For maximum cleaning effect it is recommended the cleaning solution to be preliminary heated up to 40-50 degrees centigrade.

4.2. Preparation of cleaning solution

4.2.1. Preparation of alkaline cleaning solution

Preparation of 3 % alkaline solution of Lactodaily for circulation cleaning in the milk analyser:

1. Take the package concentrated chemical Lactodaily
2. Carefully open, paying attention not to spill it.
3. In appropriate vessel (for example bucket) pour 1 l water.
4. Add the chemical and then again water up to 3 l.

Then follow the instruction for milk analyser cleaning.

4.2.2. Preparation of acidic cleaning solution

Preparation of 3 % acidic solution of Lactoweekly for circulation cleaning in the milk analyser:

1. Take the package concentrated chemical Lactoweekly
2. Carefully open, paying attention not to spill it.
3. In appropriate vessel (for example bucket) pour 1 l water.
4. Add the chemical and then again water up to 3 l.

Then follow the instruction for milk analyser cleaning.


Fig. 10 Labels for acidic cleaning chemical Lactoweekly and alkaline cleaning chemical Lactodaily



4.3. Working with the peristaltic pumps

There are two peristaltic pumps in the analyser – 1 for the sample and 1 for the cleaning.

Fig. 11 Peristaltic pump SL 16

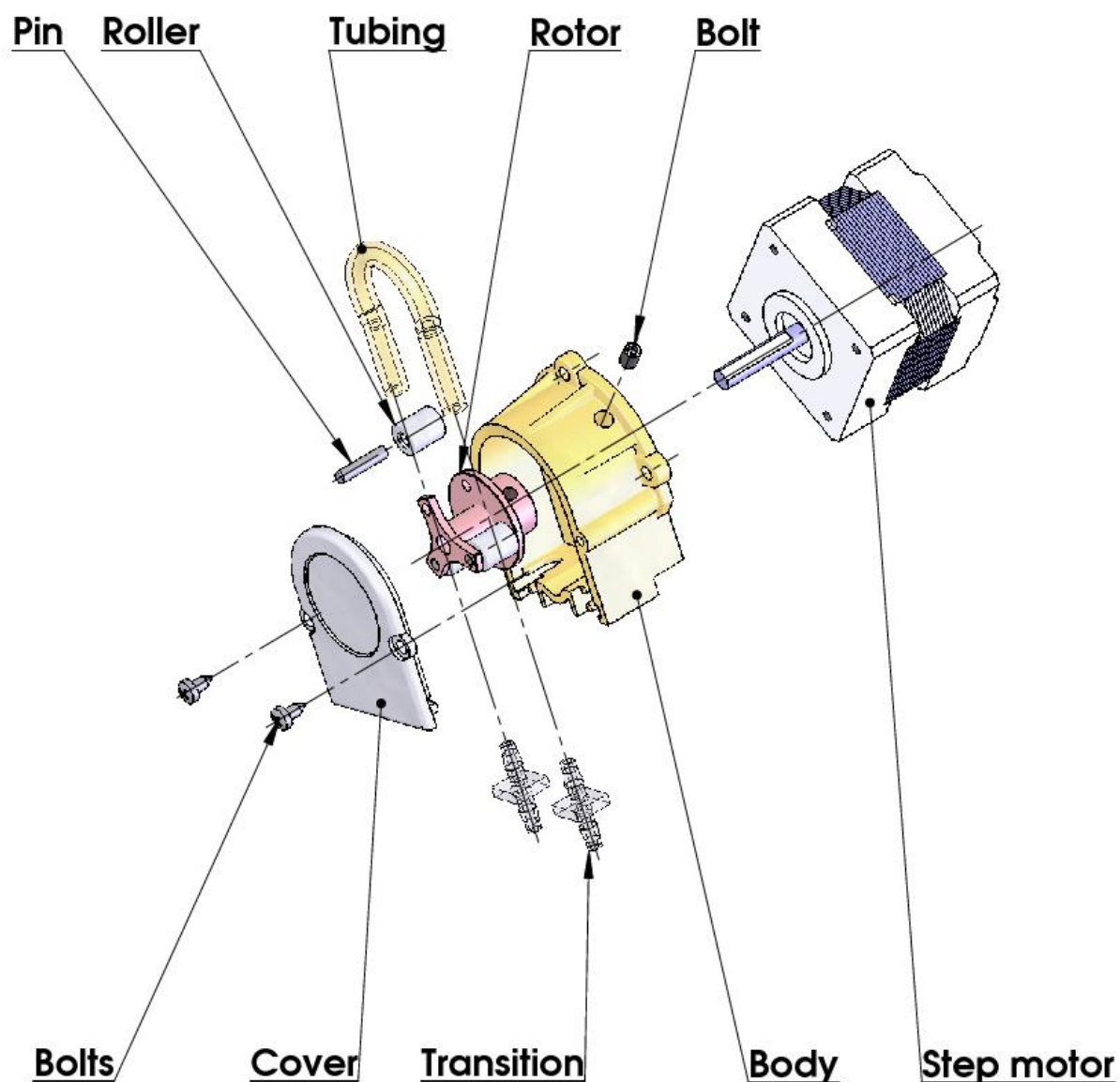

5. POSSIBLE MALFUNCTIONS AND ERROR MESSAGES, TROUBLESHOOTING

In the table below are described the possible malfunctions during the milk analyser's exploitation and ways for their repair/remedy. If the problem persists after all recommended measures are taken, please, connect the nearest service center.

Error message	Possible problem /cause	Repair/remedy
2. MA overheated Accompanied by a continuous sound signal	Overheated milk analyser	Immediately switch off the analyser. Pay attention the analyser to be situated away from direct sunlight or heating devices. Wait 5-10 minutes the device to cool down or to be normalized the ambient temperature and switch it on again.
3. Empty Camera	Insufficient quantity of the milk sample sucked in the system or air in the sample	The analyser is ready to measure the next sample. In order to avoid the future appearance of the same error message, please, check the following: <ul style="list-style-type: none"> - The sample is prepared according the instructions and there aren't air bubbles in it. - There is a real suction of the sample after starting measurement, i.e. it is obvious that the level of the milk sample in the sample holder decreases. In other case – there is damage in the suction system. - Avoid the end of the suction pipe to be above the surface of the liquid (not dipped enough). - Avoid curdling of the milk sample. Clean immediately if there is a sample curdled in the system. - In mode Measurement, after starting the measurement, remove the sample holder and see if there is no milk poured back in the sample holder.
4. Sample Overheat	Sucked overheated sample	The analyser is ready to measure the next sample. In order to avoid the future appearance of the same error message, please, check the following: <ul style="list-style-type: none"> -The sample is prepared according the instructions and its temperature does not exceed the maximum permissible sample's temperature. -Complete the procedure for checking the analyser in case of error message Empty Camera.

6. MAKING CORRECTIONS AND RECALLIBRATION OF THE DEVICE

In the process of work with the analyser there is a possibility the results to start differing between the data for some of the measuring parameters when measured with the milk analyzer and the corresponding reference method of analysis (Gerber for fat, Kjeldhal for proteins etc). In order to establish the possible discrepancy and to correct the readings of the milk analyser do the following:

6.1. Taking samples and preparation of samples for checking the accuracy of the milk analyser, making corrections and recalibration

This is a basic moment for the correct checking the accuracy of the analyser and for making correct and precise correction and calibration. It is accomplished according Appendix Sampling and preparation of samples for verification the accuracy of the milk analyzer, making corrections and recalibration.

6.2. Determination the type of the discrepancy:

6.2.1. Making measurements

Make measurements with different samples (not less than 3) with known values of a separate parameter (for example fat content), determined by the known reference methods of analysis (for example Gerber's method for determination of fat content). For more accuracy it is recommended among these samples to be also such with values, close to the lowest and highest bounds for the measured parameters.

Make 5-time measurement for each of the samples. Calculate the average value for each sample parameter, without taking into consideration the first measurement for each sample.

6.2.2. Analysing the measurement results

Make comparison between the values of the parameter from the reference sample and measured with the analyser. Make analysis of the difference received.

6.2.2.1. If the received differences are relatively constant value for samples with different content of the analysed parameter, it is necessary to make correction.

For example

M% of the reference samples:	2,20	3,00	3,80	4,60	5,20
M% average when measuring with the analyser:	<u>2,38</u>	<u>3,17</u>	<u>4,01</u>	<u>4,79</u>	<u>5,42</u>
Difference:	0,18	0,17	0,21	0,19	0,22

Conclusion: the difference is relatively constant value and correction is possible to be done with – 0,2 % (see Corrections, p6.3.3)

6.2.2.2. If the differences are not a constant value it is necessary recalibration to be done.

For example.

M% of the reference samples:	2,20	3,00	3,80	4,60	5,20
M% when measured with the analyser:	<u>2,02</u>	<u>2,93</u>	<u>3,76</u>	<u>4,75</u>	<u>5,44</u>
Difference:	-0,18	-0,07	-0,04	0,15	0,24

Conclusion: It is obvious that the difference is variable value and recalibration have to be done (See Recalibration, p.6.4).

6.3. Making corrections

6.3.1. Possible corrections, limits and changing steps

Every parameter from each calibration may be separately corrected. Below is the table with possible corrections, limits and changing steps:

Parameter	Increasing	Decreasing	Step
FAT	0.95%	0.95%	0.01%
SNF	4.75%	4.75%	0.05%
Density	4.75%	4.75%	0.05%
Lactose	0.95%	0.95%	0.01%
Salts	0.95%	0.95%	0.01%
Proteins	0.95%	0.95%	0.01%
Added water	9.00%	9.00%	1.00%
Sample's temperature	9.90°C	9.90°C	0.1°C

6.3.2. Preparing the analyzer for mode Corrections

6.3.2.1. Press the button **Enter** and without releasing it switch on the power supply of the device, wait for the starting identification messages and release the button after the following message appears on the display:

**Release button
to start setup**

After releasing the button on the display is shown:

Setup Menu

followed by possible to be entered by the operator menus:

Special modes
Corrections
Settings

Tests
pH & Co Meter
Accessories
Exit

6.3.2.2. By using buttons “up” ▲ and “down” ▼ position on **Corrections** and press **Enter**.

6.3.3. Making correction

6.3.3.1 Determining the correction mode

When starting **Corrections**, the following appears on the display:

Corrections:
Measurement
Temperature
Cond measure
Exit

Position on **Measurement** and press **Enter**. By using buttons “up”▲ and “down”▼ position on the corresponding calibration (for example **Correction 1** – cow) and press **Enter**.

6.3.3.2. Choosing correction parameter

After choosing calibration mode the display shows the following:

Cal1 Cow
Param:Fat
Correct=00.00

Edit OK Next

Using the buttons “up”▲ and “down”▼ position on the action you want to take (for example Edit) and press the button **Enter**.

6.3.3.3. Making correction

After choosing parameter (for example fat) the display shows the following:

Cal:....
Param:....
Correct= 00,00
- OK +

Using the buttons “up”▲ and “down”▼ is possible to increase or decrease the value of the measured parameter in the above pointed limits. Leaving this mode means saving the correction value and activating it.

6.3.3.4. Making verification

After the corrections are made put the milk analyser in working mode and make several times measurement of reference samples with known values of the corrected parameter. If the difference between the values of the parameter from the reference methods and milk analyser are in the limits for the parameter it may be considered that the correction is successfully made. If the discrepancy between the measurements from the milk analyser and classical methods is bigger than is necessary to make second correction according above described way.

If after the second correction the results are unsatisfactory we recommend making a calibration of the analyser. In dependence of the conditions and your requirements you may make the calibration using a personal computer

type IBM PC and the company's calibration program or autonomous - by recalibration.



When making corrections or calibrations be 100% sure in the accuracy of the reference methods result.

6.4. Recalibrating the milk analyser

6.4.1. Running the analyser in mode Recalibrate

6.4.1.1. Press the button **Enter** and without releasing it, switch on the power supply of the device, wait for the starting identification messages and release the button after the following message appears on the display:

**Release button
to start setup**

After releasing the button is displayed:

Setup menu

Followed by the possible to be entered by the operator menus:

**Special modes
Corrections
Settings
Tests**

**pH & Co Meter
Accessories
Exit**

6.4.1.2. By using buttons “up” ▲ and “down” ▼ position on **Settings** and press the button **Enter**.

6.4.1.3. In the newly opened menu, position on **Recalibrate** and press the button **Enter**.

6.4.2. Making recalibration

6.4.2.1. Choosing the calibration mode

After starting **Recalibrate**, the display shows the following:

Cal: 1 Cow		
Prev	OK	Next

You can choose the type of milk to be calibrated. By pressing ▼ (**Next**) you can switch between **Cal: 1**, **Cal: 2** or **Cal: 3**.

By pressing the button **OK** you are choosing the type of calibration.

6.4.2.2. Entering values for the separate sample parameters

The following menu is displayed:

Cal1 Cow		
FAT=f.ff		High
Edit	OK	Next

In this display the results, received by using the corresponding reference methods from *Appendix Methods* for **high-fat** milk analysis must be entered.

In this menu, with button ▼ (**Edit**) must be entered the values of the high fat milk sample

For example:

FAT=05.29

Cal1 Samp.High		
FAT=f.ff		
-	OK	+

With buttons ▼(-),▲(+) set the needed value. With next pressing of **Enter** the cursor is moved to the next number. After needed value entering completion for FAT, press **Enter (OK)** and you are going back to the previous menu:

Cal1 Cow		High
FAT=05.29		
Edit	OK	Next

With the button **▲(Next)**, choose **SNF** and in the same as above described procedure, enter the value for **SNF**. After it is finished, press **“Enter”** (OK) and you are going back to the previous menu. With button **▲ (Next)**, choose **DEN (density)** and enter the value for density; the rest of the parameters are entered in the same manner – LAC (lactose), SOL (salts), PRO (protein)

Cal1 Cow		High
PRO=f.ff		
Edit	OK	Next

If you miss to enter some of the parameters of milk, the following warning message will appear:

<p>You Must Enter Values > 00.00 Try Again</p>

Then you must press the button **Enter (OK)** and enter the missed parameters. After all the parameters are entered, press **Enter (OK)**.



You must enter values for all the measured milk parameters!!!

The screen for entering the results, received with the corresponding reference methods (See *Appendix Methods*) for the **low-fat** milk is displayed:

Cal1 Cow		Low
FAT=f.ff		
Edit	OK	Next

In the same way the values of milk with low fat sample are entered.



You must enter values for all the measured milk parameters!!!

In other case the calibration will not be correct.

6.4.2.3. Making recalibration with the available samples

After entering the values for the separate parameters of the sample, pressing **Enter (OK)** will display the following menu:

**Cal: Cow
Put sample High
5 times**

which reminds us to put 5 times the sample with high **FAT**.



The sample has to be with temperature in the boundaries 15-25°C.

Before each milk measurement stir 2-3 times the milk sample by pouring it from one vessel to another. The needed quantity is poured in the sample-holder and it is put in the recess of the analyser. Start the measurement by pressing the button **Enter**. The sample is sucked. Appears the following menu:

**Cal: Cow
Put sample: High
5 times
Temp=....**

After the sample is measured, appears the following menu:

**Cow
High
N1=..... 2=.....

Cal meas=1/5**

which reminds us to make the next measurement. Before each measurement the milk is stirred by pouring it 2-3 times from vessel to vessel. Continue the procedure till the 5th measurement.

Cal 1 Samp Low
Fat = f.ff
SNF=y.yy
Den=d.dd

Lac=l.ll
Sol=s.ss
Pro=p.pp
Exit

After 5th measurement completion automatically appears the menu, which reminds us to place the **Low fat** milk sample:

Cal: Cow
Put Sampl: Low
5 times

Stir 2-3 times the milk sample before each measurement by pouring it from one vessel to another. The needed quantity is poured in the sample-holder and it is put in the recess of the analyser. Start the measurement by pressing the button **Enter**. The sample is sucked. Appears the following menu:

Make 5 times measurement of the low FAT sample.

After 5th measurement completion automatically appears the menu:

Cal: Cow
Put sample: Water
5 times

Which reminds for 5-times water measurement.

After the 5th measurement appears the menu:

Recalibrated
Power Off-On

This means that the calibration was completed successfully and the analyzer is recalibrated for cow milk, marked as "Cal: Cow".

Switch off the power supply of the device and switch it on again.

The device is ready to work with the new calibration.

Next time when the analyser is switched on, it will be ready for work with those milk types it was just calibrated with.

If calibration with another milk type is needed, do not forget to change the calibration number for the new type of milk.

Calibration for Sheep milk will be saved as second calibration, UHT – as third. This sequence may not be followed and calibrations can be saved in whichever order is needed. Calibration can be done with different liquid dairy products using 2 representative samples.

Checking the calibration

1. Switch on the calibrated device.

2. Make sure it shows the same serial number as this already calibrated.

For checking, use the third sample with medium FAT content.

3. Measure the milk 5 times in the mode you've calibrated it.

In case that the device is not connected towards printer write down the results.

4. Ignore the first two results.

The rest three could not differ more than 0,05% FAT, 0,07% SNF, 0,7% Density one from another

7. STARTING THE DEVICE IN A SERVICE TEST/SETUP OPERATIONAL MODE. MENUS DESIGNATION

7.1. Starting the device in a service Test/Setup operational mode.

In order to start the **Setup** of the device the operator has to press the button **Enter** and without releasing it to switch on the power supply of the device, to wait for the starting identification messages and to release the button after the following message appears on the display:

**Release button
to start setup**

After releasing the button on the display is shown:

MA Setup

Followed by possible to be entered by the operator menus:

Special modes
Corrections
Settings

Tests
pH & Co Meter
Accessories
Exit

You may move in the menus by using buttons “**up**” ▲ and “**down**” ▼.

If by pressing the button **Enter** you choose a menu, each menu offers new points/submenus. When **Exit** is chosen the device leaves the **Setup** mode and returns to normal work.




Due to continuous device improvement or due to the type of ordered type device it is possible some of the functions of the analyser to be not active (options not included). In this case, when trying to enter such function will display the following message: **Not available option**

7.2. Menus Function:

7.2.1. Special modes.

Serve for choosing special (technological) working modes. After starting it the following appears on the display:

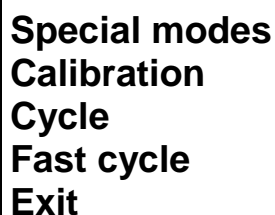

Special modes
Calibration
Cycle
Fast cycle
Exit

This mode is normally used in production conditions.

7.2.1.2. Calibration mode

In mode **Calibration** the analyser is ready to make measurement and to send the received results towards the technological milk analysers calibration system. For this purpose you need personal computer type IBM PC, company's calibration system LSC.EXE and methods for calibration of milk analysers (see the corresponding documents). To start measurement in this mode, the operator has to put a sample-holder containing milk sample in the recess of the analyser and to press the button **Enter**.

7.2.1.2. Cycle mode / Fast cycle mode

Mode **Cycle** serves for training the analysers. When you start this mode, the analyser, without additional commands, sucks the sample, makes the measurement, pours the sample out in the sample-holder and displays the received results cyclically.

7.2.2. Corrections

Serves for entering corrections in the measured data. Detailed description in point 6.3.2 and 6.3.3.


7.2.3. Settings.

Serve for assigning different working parameters (modes).

7.2.3.1. Net number.

Serves for assigning the device network number when connecting it in the production network. The possible numbers are from 0 to 15 including.

After starting this function the display shows the following:

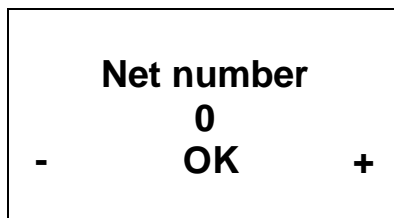

By using the button “up” ▲ the operator has the possibility to increase the number, showing the channel’s number, and by button “down” ▼, to decrease it. Pressing the button **Enter** saves the chosen channel and exits the function.



When connected in the production network each device has to have a unique number.

7.2.3.2. Recalibrate.

Serves for changing definite calibration. Methods are described in point 6.4.

7.2.3.3. Save/Rest Cal.


Through this menu you may save the new calibration in the device or to restore the old one (factory) calibration. This is necessary in case that you’ve calibrated the device for cow milk, but after that the device is not measuring correctly and you decide to restore the factory calibration settings. Position the cursor across “**Restore calibration**” and press “**Enter**”

Possibilities:


Save calibration – saves the chosen calibration in an internal buffer.

Restore calibration – restores the chosen calibration from the internal buffer.

The procedure **Save/Restore** is done for each calibration separately.

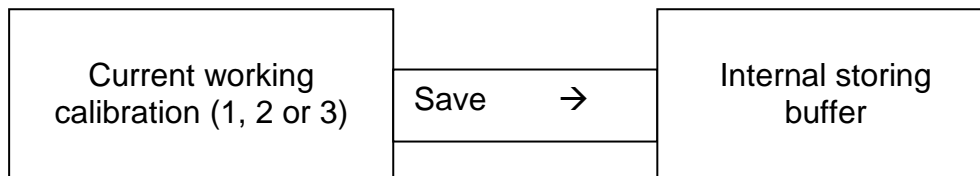

Current calibration content is not changed, the analyser continues using it, but there is a reserve copy in an internal buffer.

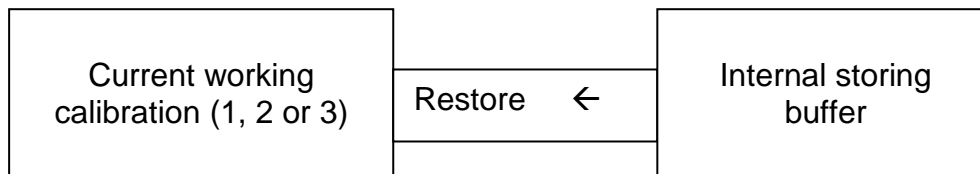


The current calibration is replaced with the calibration from the internal buffer and the analyser starts working with it. The content of the internal buffer is not changed.



If after recalibration “Save calibration” is pressed the new calibration settings will be saved over the factory settings. After that is impossible to restore the factory settings of the calibration. Save the newly made calibration only if you are sure about its correctness.

7.2.3.4. Settings Page 2.

After this menu is started the display shows the following:

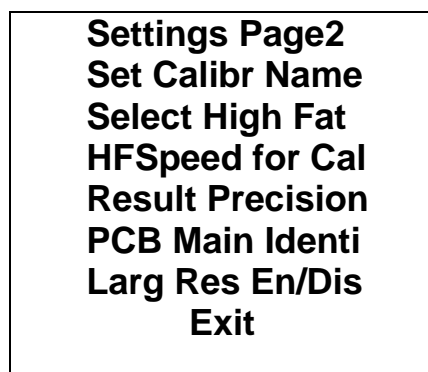

Now there is a possibility one of the following options to be set:

7.2.3.4.1. Larg Res En/Dis.

The format of the measurement data send towards the computer is set.

If the option **Large Disable**, is chosen, then only the main results are send to the computer – Fat, SNF, Density, Lac, Proteins, Added Water, sample temperature, device serial number and calibration number.

If the option **Enable**, is chosen, except the above mentioned parameters also data for Salts, Freezing Point, pH, Conductivity will be send to the computer. In this case is necessary the software in the computer to be conformable to the format of the sent data. After starting the menu, the display shows (for example):

Large Results
No

No OK Yes

7.2.3.4.2. Select High Fat

If the analyzer has embedded function for measurement of high fat products, by this menu the calibration, with which this measurement to be made is chosen. What is seen when this option is turned on is an obvious slowdown of the sample's suction speed.

7.2.3.4.3. HFSpeed for Cal

If the analyzer has embedded function for measurement of high fat products, and a new calibration for high fat measurement is needed, before starting the new calibration the operator has to start this menu. What is seen when this option is turned on is an obvious slowdown of the sample's suction speed during calibration.

Please, pay attention to the fact that switching off the power supply cancels this command action.

7.2.3.4.4. Reslt Precision

Serves for setting the precision of the measured results shown on the display. It is set separately for every parameter, the precision could be 0.01 (standard) or 0.1.

7.2.3.4.5. PCB Main Identi

Gives information about the type and the version of the analyser (LS Identity).

7.2.3.4.6. Set Calibr Name.

Sets the names of the separate calibrations. The name could be chosen from the group of predefined calibrations names or to edit a new one. When editing the new name there is a possibility all ASCII codes to be used, as letters (caps and normal), numbers and punctuation marks and popular symbols. The calibration name consists of 8 symbols.

Example:

When it is suitable to us this possibility of the analyser? For example if you have a device factory calibrated for Cow milk, Sheep Milk and UHT milk, but you need oftenly to measure camel milk. Using the methods, explained in details in Appendix Methods you may make a new calibration without need to send the analyser back to the producer for calibration. Using this procedure you may make calibrations for most oftenly analysed milk and to write down the exact calibration name, which will be shown on the display and printed on the printer.

After starting this menu the display shows:

Select Calibr
Cal1: Sheep

Exit Yes Next

There are the following possibilities:

With button **Exit** – to leave the menu.

With button **Yes** – to confirm the chosen for editing calibration name.

With button **Next** – to choose the next calibration name for editing.

If a calibration for change or edit of name is chosen, the display shows (example):

Cal:Sheep

PreDef Exit Edit

There are the following possibilities:

With button **PreDef** – to choose a calibration name from the list of preliminary given names.

With button **Exit** – to leave the menu.

With button **Edit** – to edit the new calibration name.

If a name from the preliminary given names list is chosen, the display shows:

Cal1: Sheep
UHT

Exit Yes Next

There are the following possibilities:

With button **Exit** – to leave the menu.

With button **Yes** – to confirm the chosen from the list calibration name. Now the program returns to the beginning of the menu for setting calibration names.

With button **Next** – to show the next calibration name from the list.

If it is decided a new calibration name to be edited, the display shows:

Cal1: Sheep
User Edited
Name:
Prev Set Next

There are the following possibilities:

With button **Prev** – to display the previous ASCII symbol.

With button **Set** – to confirm the ASCII symbol, shown on the display and passes to editing the next symbol from the calibration name.

With button **Next** – to show the next ASCII symbol.

After editing the last (eighth) name symbol, the display shows:

Cal1: Sheep
User Edited
Name: MilkShp
Exit Save

There are the following possibilities:

With button **Exit** – to leave the menu.

With button **Save** – to confirm already edited calibration name and to save it in the device. The program returns to the beginning of the menu for setting calibration names.

7.2.3.5. Set Base FrPnt

Through this menu you have the possibility to enter the basic freezing point separately for each calibration. For more information, see *Appendix Freezing point*. After choosing the menu the following is displayed:

```

Edit FrPoints
FrPoint Calibr1
FrPoint Calibr2
FrPoint Calibr3
-----

```

```

Exit

```

After choosing freezing point for the calibration, the following is displayed:

```

      FrPoint Calibrx
      -0.fff
-      OK      +

```

Where:

Calibrx - basic freezing point to be edited for chosen calibration;
-0.fff - basic freezing point current value;

By pressing the buttons:

"up" ▲ - you may increase the absolute value of the freezing point;
"down" ▼ - you may decrease the absolute value of the freezing point;
"Enter" - saves the edited value and exits the menu;

7.2.4. Tests.

Start different tests. Possibilities:

7.2.4.1. Test pump.

Starts pump's test. The number of the completed suction/display cycles is indicated.

7.2.4.2. Ultrasound.

Test for the ultrasonic system. Used in production conditions.

7.2.4.3. Set Amplitude.

Serves for ultrasound amplitude adjustment. It is used under production conditions or by the customer (after sensor change) according the instructions in the document SetCell.pdf.

7.2.4.4. RS232 COMPort

7.2.4.5. Keypad

7.2.4.6. USB Flash



Please, use this menu only after reading the above pointed document SetCell.pdf

7.2.5. pH meter & Co meter

Ph and conductivity measuring are additional possibilities for the analyser and are optional. Their usage is described in *Appendices PH Measurement and Conductivity Measurement*.

7.2.6. Exit

By pressing the button you may leave the program and pass towards another menu.

7.2.7. Summary of milk analyser setup functions

Analyzer Setup	
Special modes	Calibration Cycle Reports Shut Down
Corrections	Calibration 1 Fat SNF Density Lactose Salts Proteins Water Calibration 2 Fat SNF Density Lactose Salts Proteins Water Calibration 3 Fat SNF Density Lactose Salts Proteins Water
Settings	Temperature Cond measure Net number COM1 mode LCD Setup Contrast fast Contrast fine B.Light fast B.Light fine Recalibrate Calibration 1 Calibration 2 Calibration 3 Edit samp's 1 Edit samp's 2

	Edit samp's 3	
Edit FrPoints	FrPoint Calibr1	
	FrPoint Calibr2	
	FrPoint Calibr3	
Save/Rest Cal	Save Calibr 1	
	Save Calibr 2	
	Save Calibr 3	
	Rest Calibr 1	
	Rest Calibr 2	
	Rest Calibr 3	
Fan Temp		
Offs		
Settings		
Page2	Final clean cnt	
	Auto Print Res	
	Larg Res En/Dis	
	Set Calibr Name	
	Option Select	
	Sel Pump Speed	
	Speed for Cal 1	
	Speed for Cal 2	
	Speed for Cal 3	
	HFSpeed for Cal	
	OutSample Delay	
	Reslt	
	Precision	
	Fat - Precision	
	SNF - Precision	
	Den - Precision	
	Lac - Precision	
	Sol - Precision	
	Pro - Precision	
	AdW - Precision	
	Properties	
	PCB Main Identi	
	PCB LCD Identi	
Tests		
	Test pump	
	Ultrasound	
	Serial Prnter	
	KeyPad	
	Displays	
	Parallel Prn	
pH & Co Meter		

	pH	
	Calibration	
	pH	
	Measuring	
	pH	
	En/Disable	
	pH U Display	
	pH Test U[V]	
	Co Meter	
	Calibr	
	Co Meter	
	Test	
	Co Meter	
	En/Dis	
Accessories		
	Weight scale	
		Raw Test
		Calibration
		Measure
		On Line En/Dis
	Battery	
		Measure Bat U
		Ctrl Enabl/Dis
	RT Clock	
		Display Time
		Adjust Time
		Adjust Date

8. ADDITIONAL POSSIBILITIES OF THE ANALYSER

8.1. Connecting to 12 V DC power supply.

If there is a need the analyser to work on place without electrical supply available, then it could be powered by car battery or other 12 V DC external power supply. Use the 12 V power supply cable (art. number 30030).

8.2. Connecting to IBM PC

The analyser can be connected to IBM PC using the RS232 interface cable (art. number 30012, Parts and Accessories, point 16). In order to make the connection: switch off both the milk analyser and PC. Connect the RS 232 cable towards Com. Port 1 and towards the computer. Turn on both analyser and PC. Now the device is ready to communicate with IBM PC. For more detailed information regarding milk collection data program see the file resLS.

8.3. Connecting additional serial printer (option).

The interface connector for the printer is on the rear panel of the device (Com. Port 2). The printer should be connected towards it. Connect it via cables, delivered by the company-producer. If the printer is connected directly to the electrical network, then the analyser and the printer should be connected to one and the same electrical phase.

Communication parameters: 9600 bps, No parity, 8 bits, 1 stop bit. It's one-way communication (uses one line) – the analyser only sends and the printer only accepts data.

8.4. Connecting matrix printer (option).

If the device is working in milk collection center, where work load is heavy, and there are a lot of results to be printed out, we suggest an option connecting standard matrix printer (paper expenses would be less, the reliability better).

Interface connection with the printer is at the rear panel of the analyser "Parallel printer interface".

8.5. Connecting and working with external keypad (option).

It is possible external keypad (supplied by the producer of the device) to be connected to the device, using special cable to the serial printer connector.

In order to connect external keypad to the milk analyser, follow the procedure below:

1. Connect the keypad towards connector labeled Printer.
2. Connect the printer towards connector labeled **RS 232**.
3. Forward the printer towards output **RS 232**, doing the following:
 - Press the button **Enter** (on the front panel of the analyser).
 - While holding it pressed, switch on the analyser's power supply.

- Release the button **Enter**.
- Using the buttons **Up/Down** choose menu **Settings** and press **Enter**.
- Using the buttons **Up/Down** choose menu **COM1 Mode** and press **Enter**.
- With the right button (**Up**) choose **COM1 Mode – Prn**, press **Enter** and leave the menu.

With the external keypad 4 digit identification number of the milk deliverer may be entered (from 1 to 9999) and quantity of delivered milk in liters (from 0,1 to 9999.9) (accurateness up to 0.1 l). These data entering have to be completed before starting measurement. Pressing the button **Enter** on the keypad will display the following:

Enter Data
Del N:

Now enter the deliverer's identification number. After pressing **Enter** on the external keypad, the display shows:

Enter Data
Del N:xxxx
Liters =

Now, enter how many liters have been delivered and press the button **Enter** on the external keypad the following is displayed:

Del N:xxxx
Liters=yyyy.y
Are you sure?
0-No Yes-Enter

Where:

xxxx - entered by the operator identification number of the deliverer
yyyy.y - entered by the operator liters.

The operator has two possibilities:

1. To cancel the entered data by pressing the button 0 on the external keypad and to start entering them again.

2. Confirm the entered data by pressing the button **Enter** on the external keypad. The display shows currently chosen calibration. The device is ready to start measurement. After completing the measurement, the data for the deliverer are printed out.

8.6. Measuring high fat samples (option).

The standard device measures samples up to 25% fat.

On customer's request, the device could be produced with possibility to measure samples up to 50% fat. The customer can choose which calibration to have this possibility and which not, as well as during the process of exploitation to change the measuring mode i.e. to pass from measuring normal fat percentage towards high and vice versa.

What the operator sees during these passes is the difference in the speed of sucking the sample. For that purpose, the high-fat sample has to be preliminary heated up to 30C +- 3C.

To choose the mode, follow the sequence below:

Setup->Settings->Settings

Page2->Option

Select->SelfPumpSpeed->Speed for Cal x

After which the display shows:

Calibr x Pump Speed
Normal OK HiFat

By pressing the correspondent buttons the operator can choose the type of measurement and to exit the menu.

When changing the type of measurement on a calibration is necessary a new calibration of the device on the new speed to be done. When calibrating measuring high fat sample, before starting the calibration procedure, the operator has to choose from the menu:

Setup->Settings->Settings

Page2->Option

Select->SelfPumpSpeed-> HFSpeed for Cal

By which the device passes in a mode of measuring high fat samples. This calibration mode is active till the power supply of the device is switched off i.e. it has to be always set if the device will be calibrated for high fat measurement.

8.7. Working with embedded accumulator (option).

On customers' request accumulator could be embedded in the milk analyser. In this way it could work independently on the electrical network. In this case the Power switch on the rear panel of the analyser serves for close / open the electric circuit of the outer power supply. After closing the circuit, in order to start the analyser, the operator has to press the Start button on the front panel of the analyser.

Switching off the analyser could be done on one of the following ways:
Through the analyser – by switching the Power button on the rear panel of the analyser.

Through the software – by using the option Shut Down from the Calibration and Working mode menus.

After pressing and holding pressed the button Enter (when choosing calibration), the following is displayed:

Milk selector
Cal1 – Cow
Cal2 – Sheep
Cal3 – UHT

Cleaning
Final clean
Shut Down

By choosing the option Shut Down the operator may switch off the power supply of the analysers.

In case that the external power supply from the electrical network, the analyzer automatically starts working with power supply from the battery. Possible working modes with battery:

1. The analyser is powered by the battery, works normally, the display shows the following message

Battery Powered

2. The analyser is in mode Sleep – do not accept any commands, expects change in the power (supply from switching adapter for starting the battery), display shows the following

Low Battery

3. The battery power is under the allowed minimum, milk analyser is automatically switched off. Before switching off, for a short time, the display shows


Out of Battery

8.8. Embedded real time clock (option).

On customer's request, a real time clock could be embedded in the device, showing astronomical time and date. The clock is powered by battery, so it is independent on power supply of the device. The advantage of this option is that on the print out with the measurement results are shown also the exact time, when the measurement is made, for example:

```
Time: 18:14:33
Date: 22:03:2007
Lactoscan MCC30
Serial Number:0002
Calibr 1 Camel
Results:
Fat.....00.00%
SNF.....00.00%
Lactose.....00.00%
Solids.....00.00%
Protein.....00.00%
Temp. Sample...21.0°C
```

The embedded clock is controlled by the device's Setup, from the main menu Accessories, submenu RT Clock. When chosen, the display shows the following:



Display Time
Adjust Time
Adjust Date
Exit

Using these menus, the operator has the possibility to show on the display the current time and date, and if necessary – to adjust them.

The time is shown in the format:

hh:mm:ss
dd:mm:yy

Exit

where:

- hh - the current hour
- mm - current minutes
- ss - current seconds
- dd - current day
- mm - current month
- yy - current year

9.11. USB Flash-drive (pen drive) (option).

This option allows connecting standard USB Flash Drive (Pen Drive) towards the Lactoscan. Its usage allows automatically (without operator's interfere) creating total archive of all measurements, made with the analyzer. Actions necessary for switching on/off the USB Flash Drive are analogous to the same when working with PC. After USB Flash Drive is inserted in the coupling at the front panel of the analyzer, Lactoscan starts writing there the results from each measurement. If you do not need to make an archive, then simply remove the flash drive from the analyser. As when working in PC, when the flash reads/writes its LED blinks. In this case wait for the LED stop blinking and then remove the flash drive from the analyser.

Data, collected from measurements, are saved in MS Excel file. The name of the file is formed by the date, when the measurement is made (in case there's embedded clock option). For example – 111213.xls means, that the file was created in year 2011 (last two digits from the year), month 12 –

December, day13. Tabs with results are in accordance with software resLS, i.e. the file, created on the flash drive is with the same structure and may be processed with the same, already existent client software, used for collecting data with PC (see the manual of the software resLS.exe). Meaning of results tabs are in the file LS-USB.XLS.

The difference is, that if there's flash drive connected to Lactoscan during measurements, there's no need the analyzer to be connected to PC all the time. In the presence of more than 1 analyzer, each of them may be equipped with flash drive and at the end of the shift (working day or after a definite period of time (chosen by the customer) files from different analyzers to be processed on only one PC.



It is recommended this option to be used in combination with RTC and Keypad. Availability of clock option ensures archive of the time when the measurement was taken. Clock adjustment can be made with the software tool LS-AdjustRTClock.exe (see its manual). In presence of Keypad, the deliverer's number can be entered (or some kind of number/identifier of the sample) – number from 1 to 9999 and liters/kg of the milk.

There's no requirement for type and producer of the Flash Drive, as well as its storage capacity. Other files, working under MS Windows may be also present on the flash drive.

Approximately, results from one measurement, take less than 100 bytes. For example, flash drive of 8 GB, can store more than 80 000 000 measurements, which makes it unlimited from point of view archive of measurement results.



Use this option ONLY with flash-drives. Do not connect any other USB devices.

APPENDICES**APPENDIX 1: PREPARATION OF SAMPLES FOR MILKANALYSERS' CALIBRATION**

For calibration are needed samples of cow milk with the following parameters:

		Low Fat	High Fat	Middle
1	Cow	2,2%	5,2%	3,6%

For the calibration are needed:

1. Distilled water
2. Min. 3 milk samples with known values for fat, SNF, protein, density, lactose, salts.

Calibration samples have to be with low, middle and high values of the analyzed components. Samples have to be representative for given milk type. Volume of the sample has to be enough for making min 5 measurements for each sample – not less than 1,00 l. Changes in the analyzed parameters in the samples have, if possible, to cover the whole measuring range – i.e. used samples to be with low, middle and high content of the analyzed components.

Methods of milk samples preparation for calibration.

For milk sample with middle value of the analysed components we recommend to use milk taken from not less than 10 animals from most common in the region breed.

Sample with low and high value are prepared on the following way:

1. Pour the fresh milk with FAT at about 3.7% in a separating funnel.
2. Leave the funnel with the milk in refrigerator for 12 hours at temperature +5-+8 ° C.
3. Draw the substratum of the separated milk in a vessel, mix it well, pour it and heat it in water-bath up to 20°C.
4. Pour the upper layer in another vessel.
5. Determine the concentration of the measured components (FAT, protein, SNF, density, lactose, solids) by using certified methods.



The analyser's accuracy depends only on the correctness of the chemical analysis of the components in the samples and the normal acidity during calibration!

It is recommended the first cow milk sample with low fat content to be with the following parameters:

2-2,3% FAT; 8.7-9% SNF; 3,3-3,5 % Protein; 4,8-4,9% Lactose; 0,75 Salts; 1030-1033 kg/m³ Density.

The second cow milk sample with high fat content to be with the following parameters:

5-5,3% FAT; 8.4-8,79% SNF; 3,1-3,2% Protein; 4,6-4,7% Lactose; 0,7 Salts; 1028-1029 kg/m³ Density.

If, after milk's separation you do not obtain samples in the requested range, then, by adding milk with high fat value into the low fat milk sample you can obtain necessary value-2,3%

Analogous to this, by adding low fat milk sample into a milk sample with high fat value you may receive 5,3%

Samples with medium values are received by mixing low fat and high fat samples in necessary proportion.

If there is a need of longer sample storing they have to be preserved; the most commonly used preservative is potassium dichromate (K₂Cr₂O₇) - 1 g for 1 000 ml.

When using samples, stored shortly, preliminary pour the sample from one vessel to another in order to distribute the milk components evenly paying attention not to form foam in the sample.

When the samples are stored for a longer period it is recommended to warm it up to 35-45 °C, and the vessel to be shaken carefully. In case that there is a cream stuck on the vessel's surfaces – remove it. The sample is poured from vessel to vessel several times and is cooled down (advisable to 20 °C /.



If there is separated liquefied fat or white particles with irregular form on the vessel's walls reliable results could not be received.

Because it is very difficult both lactose and salts to be measured but they are substantial and influence in great extend when determine added water. That's why it is better both lactose and salts to be calculated by using SNF results. The milk must be for sure without added water.

If you are unable to make the analysis of milk in certified methods in a pinch you can use the following formulas:



DETERMINATION OF THE BASIC PARAMETERS IN THE MILK SAMPLE BY USING FORMULAS IS NOT AS PRECISE AS USING THE ARBITRARY METHODS, BUT IS SUITABLE FOR USAGE IN FIELD WORK.

1. Determination some of the parameters by formulas

There is dependence between the different parameters in milk and its density, which may be expressed with mathematical equation. On this base different formula, tested and confirmed by the classical laboratory methods for analysis, are developed. We recommend the following:

2. SNF determination.

For determination of SNF the correlation dependence exists between the milk's density, fat and SNF in the milk. When the density and the fat are known, the SNF can be calculated.

There are several formulas with different applicability.

A/ When the salts and fat are known

SNF is calculated by subtracting the fat percentage from the salts.

$$SNF = \text{Salts} - F (\%)$$

Where

Salts – salts in (%),

F – fat content in (%),

This formula is used for determination of SNF in whey, buttermilk, and cream.

B/ Known quantity of fat and density (most commonly used method when maximum accuracy is needed).

We recommend the following formula:

$$SNF = \frac{0,075 * F\% + 100 - 100 / \text{density}}{0,378}$$

This is a universal formula and actual for milk of almost all kind of cows and sheep all over the world.

3. Determination of lactose content

We recommend the following formulas:

A/ for cow milk

$$\text{Lact.} = SNF * 0,55 (\%)$$

Where

SNF – content of SNF in percentages (%),

0,55 – constant coefficient.

B/ for sheep milk

$$\text{Lact.} = SNF * 0,45 (\%)$$

Where

SNF – solids-non-fat content in percentages (%),
0,45 – constant coefficient.



This is an actual coefficient for sheep breeds on the territory of the Balkan Peninsula.

4. Determination of salts content

We recommend using the following formulas:

A/ for cow milk

$$\text{Salts} = \text{SNF} * 0,083 (\%)$$

Where

SNF – solids-non-fat content in percentages (%),
0,083 – constant coefficient.

B/ for sheep milk

$$\text{Salts} = \text{SNF} * 0,075 (\%)$$

Where

SNF – solids-non-fat content in percentages (%),
0,075 – constant coefficient.



This is an actual coefficient for sheep breeds on the territory of the Balkan Peninsula.

5. Determination of total proteins content

We recommend using the following formulas:

A/ for cow milk

$$\text{Protein} = \text{SNF} * 0,367 (\%)$$

Where

SNF – solids-non-fat content in percentages (%),
0,367 – constant coefficient.

B/ for sheep milk

$$\text{Protein} = \text{SNF} * 0,475 (\%)$$

Where

SNF – solids-non-fat content in percentages (%),

0,475 - constant coefficient.



This is an actual coefficient for sheep breeds on the territory of the Balkan Peninsula.

APPENDIX 2 FREEZING POINT DETERMINATION

1. Methods for determination.

The milk analyser determines the freezing point of each sample and the quantity of added water. The milk analyser does not measure the freezing point, but calculates it from the components it depends on. The basic components in the milk are water, solids, lactose, FAT, proteins, minerals (salts) and acids. The freezing point depends only on the diluted in the milk components and quantity of the solvent (in the milk it is water). The ultrasonic technology allows direct measurement of FAT, proteins, lactose + salts (the soluble components, only influencing the freezing point), and the quantity of the solvent in % is determined by $100\% - \text{total solids \%}$, $\text{total solids} = \text{lactose \%} + \text{FAT \%} + \text{proteins \%} + \text{salts \%} + \text{acids \%}$.

Without understanding the meaning of the freezing point – determined or shown from the milk analyser added water result easily may lead to a mistake for the value of this parameter.

2. The basic freezing point.

Milk freezes at lower temperature than water. The average freezing point of the raw milk in the most regions is at about $-0,540^{\circ}\text{C}$. The average reading for your region is called “basic” freezing point.

The freezing point of milk is a “physiological constant”. This does not mean that it will not vary. In fact feed, breed, season, time of lactation, climate, whether the sample is taken at the beginning, middle or end of lactation – all these factors will have an effect on the freezing point of the individual sample. This means that there is an average value of all these numbers. The more samples used in obtaining this average, the more reliable it is as a base. Or the basic freezing point is an average of freezing points of milk, taken from many cows. When a laboratory checks a producer, it is only comparing the average of the producer’s cows against a larger area average.

The Health authorities establish the basic freezing point or agriculture departments in some regions, sometimes by universities, separate dairy producers, or their associations. Frequently, tolerances have been established on top of a basic freezing point to allow some variations in the milk as well as device or operator variations.

Without mentioning the basic freezing point, the Association of Official Analytical Chemists now recommends an upper limit freezing point at $-0,525^{\circ}\text{C}$ (2,326 standard deviations above the most recently determined North American average of $-0,5404^{\circ}\text{C}$), below which there will be at 95%

confidence that will show 99% of all freezing point determinations on milk without added water:

“if the freezing point is $-0,525^{\circ}\text{C}$ or below, milk may be presumed to be free of water or may be confirmed as water free by tests, specified below. If the freezing point is above $-0,525^{\circ}\text{C}$, milk will be designated as “presumptive added water” and will be confirmed as added water or added water free by tests specified below. Evaluate extreme daily fluctuations in the freezing point of herd, pooled herd, or processed milk for presence of added water”.

“Presumed added water”, as described above, must be “confirmed” by means of tests on authentic milk samples obtained as specified in the AOAC METHODS.

After determination the freezing point of your sample via the milk analyser, the added water is calculated using the following formula:

$$\text{AddedWater} = \frac{\text{FrPoint}_{\text{Base}} - \text{FrPoint}_{\text{Calc}}}{\text{FrPoint}_{\text{Base}}} * 100[\%]$$

Where:

FrPointBase is the basic freezing point

FrPointCalc is measured freezing point

Note:

If the freezing point is not correctly determined, the result for the added water is not valid. In this case results for FrPoint and AddWater are not shown on the display and on the printout from the printer. If the density of the measured sample is 0, the result for AddWater is not valid and is also not shown on the display and the printouts.

Sample:

First variant

If you've entered for milk analyser basic freezing point -0.520°C (according article 5.9 of the EU Milk Hygiene Directive 92/46/EEC), measured freezing point -0.540°C , using the above pointed formula you'll receive $-3,8\%$. Because it is not possible the added water to be negative value, the milk analyser indicates 0% added water. The reason for this is the tolerance in the basic freezing point, reasons for which are described below.

If in the same milk we add 3,8% water, and the basic freezing point is the same, the milk analyser will measure freezing point -0.520°C , and will indicate again 0% added water.

Second variant

If you've entered for the device basic freezing point -0.540°C , measured freezing point -0.540°C , the milk analyser will indicate 0%. When you add 3,8% water, the device will indicate 3,8%-added water.

From the above mentioned follows that it is very important to enter correct basic freezing point in the device.



The device's results for added water may give information about doubt of added water in the milk and the exact value of this added water may be determined after a "cowshed sample" is taken and the result for the freezing point, measured by the milk analyser of the "cowshed sample" is entered as basic freezing point in the formula for calculation of added water.

Then the result from this formula will give us the absolute value of the added water for the corresponding milk supplier.

APPENDIX 3 PH MEASURING (ONLY FOR DEVICES WITH EMBEDDED FUNCTION)

1. General information

PH probe is a unit, measuring the solution acidity or alkalinity degree. It is measured on scale of 0 to 14. The term pH is derived from "p", the mathematical symbol for the negative logarithm, and "H", the chemical symbol of Hydrogen. The formal definition of pH is the negative logarithm of the Hydrogen ion activity.

2. pH Electrode

For pH measurement the milk analyser needs a combination electrode, compatible with most pH electrodes that have BNC connectors and zero potential (the pH where the mill volt output of the electrode equals 0) near 7 pH.

2.1. Electrode part

The electrode is the most important part of the pH measurement. The electrode glass membrane is fragile and must be handled with care. To protect the glass membrane and to maintain activation, a protective rubber cap containing a suitable storage solution covers the glass membrane.

2.2. Electrode care & Electrode maintenance

pH Electrodes are susceptible to dirt and contamination and need to be clean regularly depending on the extent and condition of use. At no time should one touch or rub the glass bulb as this causes the build-up of electrostatic charge.

2.3. Storage

For best results, always keep the pH bulb wet. An optimal storage solution for combination electrode is pH 4 buffer with 225 grams of KCl per liter. Table salt, NaCl, can be used if KCl is not really available. Other pH buffers or tap water are also acceptable storage media, but avoid storage in de-ionized water. The protective rubber cap filled with the buffer solution provides ideal storage for long periods.

2.4. After Use

After measurement is completed, follow the sequence below for storage.

- Wash the electrode and reference junction in de-ionized water.
- Close the refilling hole by returning its rubber sleeve or stopper cap. (Necessary for only refillable electrode).

- Store the electrode as mentioned above (see section Storage).

2.5. Electrolyte Replacement (for refillable electrode only).

The reference electrolyte needs to be refilled when the electrode has been used for a long period, or when the internal electrolyte has dried up. To accomplish this, follow the procedure described below.

- Remove the protective rubber cap or sleeve;
- Remove the protective rubber sleeve to expose the filling port of the electrode;
- Remove the old reference electrolyte with a syringe;
- Fill the new reference electrolyte.

2.6. New electrolyte preparation:

- Open the KCl container;
- Add in de-ionized water until it reaches the level of 20 ml;
- Close the container and shake it to dissolve the KCl;
- Add in fresh electrolyte until it reaches the level of the refilling port. The reference electrolyte used should be 3M(Mol) KCl;
- Replace the rubber sleeve.

2.7. Re-use the electrode.

- Rinse the liquid junction with de-ionized water.



If these steps fail to restore normal electrode response, you may attempt to rejuvenate it (See: Electrode Rejuvenation).

2.8. Electrode cleaning

Electrodes which are mechanically intact can often be restored to normal performance by one or combination of the following procedures.

- Salt deposits:

Dissolve the deposit by immersing the electrode in tap water for ten to fifteen minutes. Then thoroughly rinse with de-ionized water. Wash the electrode pH bulb in a little detergent and water. Rinse electrode tip in with de-ionized water.

- Oil/Grease films:

Wash electrode pH bulb in a little detergent and water. Rinse electrode tip with de-ionized water.

- Clogged Reference Junction:

pH electrodes have junction, which allows the internal fill solution of the measuring electrode to leak out into the solution being measured. The

junction can become clogged by contamination in the solution. If a clogged junction is suspected it is best to clear the junction.

Heat up the diluted KCl solution to 60-80°C. Place the sensing part of the pH electrode into the heated KCl solution for approximately 10 minutes. Allow the electrode to cool while immersed in some unheated KCl solution.

- Protein Deposits

Prepare 1% pepsin solution in 0.1 M HCl. Allow the electrode to stand in this solution for five to ten minutes. Rinse the electrode with de-ionized water.

2.9. Electrode activation

Generally, if the procedure of storage and maintenance had been closely followed, the electrode can

be used immediately. However, should the electrode response become sluggish, it may be possible that the bulb has dehydrated.

The bulb can be dehydrated by immersing the electrode in an ideal storage solution (e.g. buffer pH 4 solution) for 1-2 hours. If this fails, the electrode may require re-activation. If the above procedure does not reactivate the electrode to acceptable status, try rejuvenation the electrode by following the procedure outlined below.

2.10. Rejuvenation Procedure

Dip and stir the electrode in freon or alcohol for 5 minutes.

Leave the electrode in tap water for 15 minutes.

Dip and stir the electrode in concentrated acid (HCl, H₂S₄) for 5 minutes.

Leave the electrode in tap water for 15 minutes.

Dip and stir in strong base (NaOH) for 5 minutes.

Leave the electrode in tap water for 15 minutes.

Test with standard calibration solution.

Finally, test with standard calibration buffer solution to see if the electrode yields acceptable results. You may repeat again for better response (maximum 3 times). If the response does not improve, then the electrode has completed its useful life. Replace with a new electrode.

2.11. Electrode Lifespan

pH electrodes have a finite lifespan due to their inherent properties. How long a pH electrode will last will depend on how it is cared and the solution it is used to measure. Even if an electrode is not used it still ages. Electrode demise can usually be characterized by a sluggish response, erratic readings or a reading, which will not change. When this occurs an electrode can no longer be calibrated. pH electrodes are fragile and have a limited lifespan.

How long an electrode will last is determined by how well it is maintained and the pH application. The harsher the system, the shorter the lifespan. For this reason it is always a good idea to have a back-up electrode on hand to avoid any system down time.

3. Buffer Solutions

Buffers are solutions that have constant pH values and the ability to resist changes in that pH level. They are used to calibrate pH measurement system.

PH buffer solution description (Pharmacopoeia standard)

Use only these types' standard buffers for calibration!

Description	pH 7.00±0,01/20°C	pH 4.00±0,01/20°C
Composition	Potassium dihydrogen phosphate, Di-sodium hydrogen phosphate	Borax, Sodium hydroxide solution
Temperature parameters	10°C - 7.06 25°C - 6.99 20°C - 7.00 30°C - 6.98 40°C - 6.95 50°C - 6.91	10°C - 4.00 25°C - 4.00 20°C - 4.00 30°C - 4.00 40°C - 4.00 50°C - 4.05

4. pH Electrode Calibration

pH Electrodes are like batteries; they run down with time and use. As an electrode ages, its glass changes resistance. For this reason, electrodes need to be calibrated on a regular basis. Calibration in pH buffer solution corrects for this change.

Calibration is an important part of electrode maintenance. This assures not only that the electrode is behaving properly but that the system is operating correctly.

Usually pH meters require calibration at 2 specific pH values. One calibration is usually performed at pH 7, second pH 5.



It is the best to select a buffer as close as possible to the actual pH value of the sample to be measured. Use standard calibration buffers that the temperature and the sample solution are the same.
Use the operation manual for the corresponding pH meter.

For pH electrodes, originally supplied with milk analyser read the following information:

Temperature compensations

The output of pH electrodes varies with temperature in manner, predicted by theory. When needed, Sensorex can supply electrode holders with build-in automatic temperature compensators. The need of automatic compensation depends on the temperature variation, the pH value being measured. At pH of about 7 there is no error due to temperature and, of course, at a constant temperature there is no error. As shown in the following table, the pH error due to temperature is a function of both the temperature and the pH value being measured. At a pH of about 7 there is no error due to temperature and, of course, at a constant temperature there is no error. The more the temperature changes from the ambient calibration temperature and the more the pH departs from 7 the greater is the pH error.

pH temperature error table

°C	pH										
	2	3	4	5	6	7	8	9	10	11	12
5	.30	.24	.18	.12	.06	0	.06	.12	.18	.24	.30
15	.15	.12	.09	.06	.03	0	.03	.06	.09	.12	.15
25	0	0	0	0	0	0	0	0	0		0
35	.15	.12	.09	.06	.03	0	.03	.06	.09	.12	.15
45	.30	.24	.18	.12	.06	0	.06	.12	.18	.24	.30
55	.45	.36	.27	.18	.09	0	.09	.18	.27	.36	.45
65	.60	.48	.36	.24	.12	0	.12	.24	.36	.48	.60
75	.75	.60	.45	.30	.15	0	.15	.30	.45	.60	.75
85	.90	.72	.54	.36	.18	0	.18	.36	.54	.72	.90

0 pH Error Range

Less than .1 pH Error Range

5. PH helpful hints

For greatest accuracy in pH measurement, follow these guidelines:

- Use the same technique to measure samples, which was used for calibration.
- Be consistent with stirring rates, times and conditions.
- Calibrate with buffers, which are close in temperature to that of the sample.
- Calibrate the pH electrode regularly, e.g. once a day for accuracy to within 0.05 pH.

- Use fresh buffers for calibrations. Avoid contamination of the stock buffer solution and do not use it beyond the expiry date.
- Keep all connections dry.
- Immerse the electrode far enough into the solution to insure the reference junction is below the surface.
- Allow adequate time for the electrode to stabilize in standards and samples before taking a reading.
- Clean the electrode periodically. Allow more time for aged electrodes.
- Do not use the pH electrode in solutions of fluoride ion at low pH. This will etch the glass membrane.
- Sulphide vapors can permeate the electrode wick and contaminate the reference element. Minimize contact in such environments and change the reference electrolyte frequently.

Milk acidity and pH

In the table below the connection between pH value and acidity in °T is shown

°T	deviation	pH average value	°T	deviation	pH average value
Raw milk			Pasteurized milk		
16	6,74-6,70	6,72	16	6,68-6,64	6,66
17	6,69-6,65	6,68	17	6,63-6,58	6,61
18	6,64-6,58	6,62	18	6,57-6,52	6,55
19	6,57-6,52	6,55	19	6,51-6,46	6,49
20	6,51-6,46	6,49	20	6,45-6,40	6,43
21	6,45-6,40	6,43	21	6,39-6,35	6,37
22	6,39-6,35	6,37	22	6,34-6,30	6,32
23	6,34-6,30	6,32	23	6,29-6,24	6,26
24	6,29-6,24	6,25	24	6,23-6,19	6,21

1. Preparation for pH measurement

When the analyzer is with pH measuring option, it is received from the customer with pH probe packed separately and there's a stopper on its place. If you need to measure pH follow the procedure below:

1. Loosen the nut anti-clockwise.
2. pull up the stopper

Fig. 11 Placing the stopper




3. Carefully place the pH probe paying attention not to remove the sealing O-ring.

1. Place the probe with the nut in the hole and tighten it.



It is very important to close the nut tightly, paying attention not to allow air to enter the system.

Fig. 12 Placing the probe

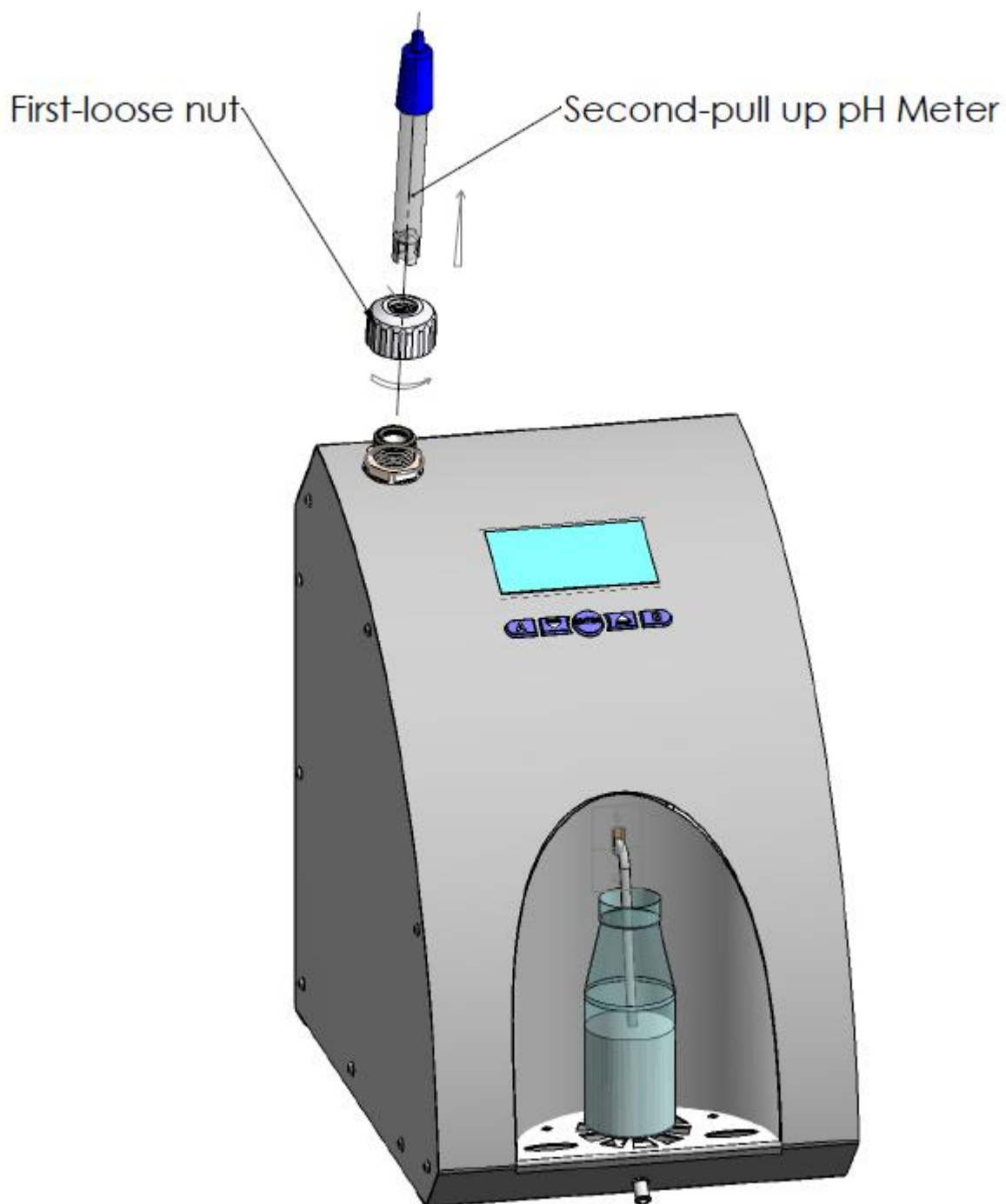




If you are working with the analyser regularly (each day) do not remove the probe after work..



If you 'll not use the analyser more than 2 days, you must take out the probe and to place the stopper back.

The pH probe must be stored separately as per the instructions of point 2.3 Appendix 3.

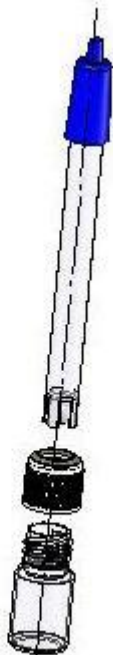

7. PH measuring.

Measuring pH is an additional feature of the analyser and is optional.

Remove the protective rubber cap of the pH electrode. Take care to handle it appropriate in order not to be damaged. Use de-ionized or distilled water to rinse the electrode before usage. Fill in the sample holder with milk, put it in the recess of the analyser and dip the pH electrode into the milk sample, ensuring complete dip of the electrode in the sample. Stir gently for homogenization of the sample.

Measuring can be done in two modes:

Off line by starting the menu **pH & Co Meter | Measuring**, when the analyser works only as a pH meter.

On line automatic pH measuring, when measuring the rest of the sample's parameters.



When starting work with pH meter first connect the probe/sensor, and then the power supply of the device.

Having in mind the characteristics of the process of pH measuring it is necessary to dip the pH probe in the sample and then to press the button Enter.

pH probe is knee-joint (mobile) and is adjusted towards the input pipe of the analyser.


Move them forward in order to be possible to dip probes' end in the milk sample and place the sample holder in the recess of the analyser.

After starting the menu **pH & Co Meter** the following message appears on the display:

```
pH Calibration
pH Measuring
pH En/Disable
pH U Display
-----
pH Test
Co Meter Calibr
Co Meter Test
Co Meter En/Dis
-----
Exit
```

8. pH Calibration.

Fig. 13 Calibration buffers

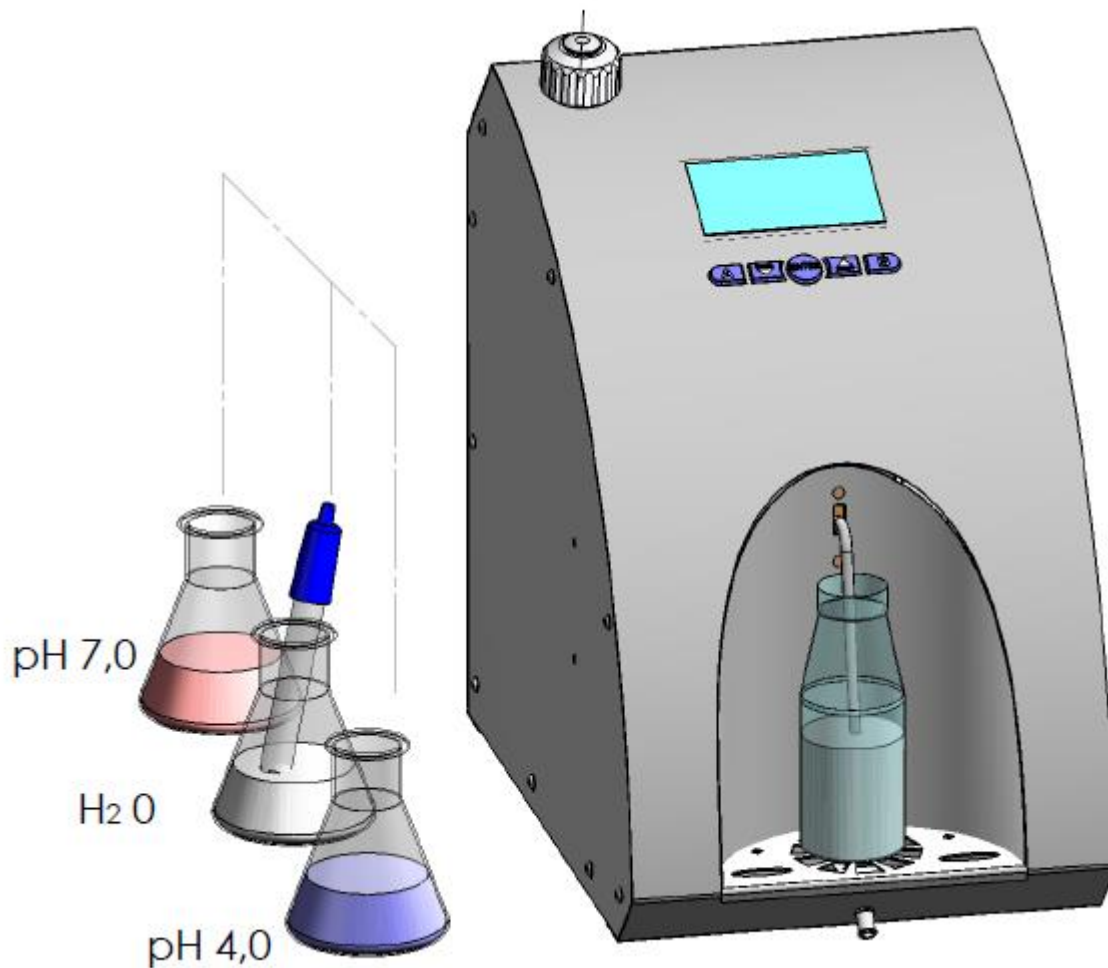

Take out the sensor. You may not place the nut, because the analyzer do not sucks during calibration.

Serves for calibrating the pH meter. For that purpose use 2 standard buffers, shown on the display as Low buffer (for example 5.00 pH) and High buffer (for example 7.00 pH). The procedure is following:

1. Start menu **Calibration**.
2. Put the pH probe in Low buffer.
3. By pressing the buttons of the analyser enter the exact buffer value. The display shows:

pH Calibr
Put Izopot buff
Buf=xx.xxx

4. Input the buffer value when the probe is in its isopotential point and with which our probes are accompanied with.

Then the display shows:

pH Calibr
Put Izopot buff
Buf=xx.xxx
V=x.xxxV Set

Where **x.xxx** is measured by the probe voltage.

5. After the readings calm down, the operator has to press the button below the inscription **Set**.

6. Repeat the procedure with **Next buffer**.

The following is displayed:

pH Calibr OK

This means that the procedure for analyser calibration was successfully completed.

7. pH measuring mode is started automatically.

8. Check calibration correctness by measuring buffer solution 7.00.

9. pH Measuring.

After starting this menu the measurement is done in mode off line, i.e. the analyser works only as a pH meter. The operator has to dip the probe in the sample and on the display the following is shown:

pH measuring
x.xxxV
y.yy pH
Exit

Where:

x.xxx – measured by the probe voltage

y.yy – measured probe's pH

By pressing the button **Exit**, the operator may exit the program and to pass towards another menu.

10. pH En/Disable

Serves for enabling/disabling the pH measuring during normal work of the analyser - On line. After starting it the display shows:

pH Measuring		
XXX		
No	OK	Yes

Where **XXX** is the current situation of the working mode. By pressing the buttons below the corresponding inscriptions it could be changed, as **Yes** – means that during normal work of the analyser – measuring the rest of the parameters, pH will also be measured. If **No** is chosen, then pH is not measured.

11. pH U Display

Serves for allowing/forbidding the value of the pH probe voltage during pH parameter measuring. After starting it the display shows:

PHUDisplay		
XXX		
No	OK	Yes

Where **XXX** is the current state of the displaying mode. By pressing the buttons below the inscriptions it could be changed, as **Yes** – means that during pH measuring the voltage of the pH probe will be shown. If **No** is chosen, it will not be shown. It refers to both of the measuring modes.

12. pH test

serves for testing the measuring system in production mode.

APPENDIX 4 CONDUCTIVITY MEASURING (ONLY FOR DEVICES WITH EMBEDDED FUNCTION)

1. Method of determination.

Conductivity (or Electrolytic Conductivity) is defined as the ability of a substance to conduct electrical current. It is the reciprocal of the resistance.

In a healthy animal*, the mean value of electric conductivity is:

Milk type	Conductivity values
Cow milk	between 4 to 6 mS/cm (18°C);
Sheep milk	between 3 to 5 mS/cm (18°C);
Buffalo	between 2,5 to 5 mS/cm (18°C);

*These values depend on the geographical region, the breed and on other factors.

Milk conductivity changes on the concentration of ions in the milk:

Added water, sugar, proteins, insoluble solids	Decrease the ion's concentration. Milk conductivity decreases.
Added salts	Increase the ion's concentration. Milk conductivity increases. Often the milk is falsified by adding salt: towards milk with good characteristics: fat 4%, SNF 8,8, conductivity 4,5 are added salt and water. Then the results are changed to 3,2 and 8,8, conductivity 10. In other words adding water regulates the increased value of SNF and density till normal (within the boundaries/parameters) and even the fat is normal. By the values of these parameters may be determined if the sample is falsified, but the only characteristic, proving this is conductivity, which is out of boundaries nevertheless added water. But be careful, as the falsification is not the only possible reason for conductivity increasing. The other possibility is mastitis that's why we recommend using another (chemical) method for checking it.
Significantly extreme value (6,5 - 13,00 mS/cm (18°C))	Should indicate the development of mastitis. Infections damage the tissue of the udder. This allows sodium and chlorine ions from the blood to

	be released into the milk. The concentration of ions in the milk is thereby raised, and it can more easily conduct an electrical current - the conductivity of the milk increases.
--	------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

Milk conductivity can be used as tests for degree of water evaporation in condense milk production.

Milk conductivity change notifies of powder (dry) milk solution rate.

2. Conductivity measurement

Conductivity measurement is additional possibility of the analyser and is delivered on customers request/

3. Co Meter Calibr

Serves for conductivity measuring system calibration. Clean the analyzer before starting conductivity measurement. (see p. 4.1). You need a standard buffer with conductivity 5.02[mS/cm] (you may order it for delivery together with the analyser), with temperature 18°C. After starting this mode, the analyzer makes preparation for measurement and when it is ready, the following message is displayed:

Co Meter Value Base= 5.02

The basic value of the buffer solution is shown. We can use buffer solution with another conductivity (from 4 to 5 mS/cm) or refer to the Note at the end of this point. The value of the used buffer must be changed or confirmed. Pressing ENTER confirms, + or – increase or decrease the value. After 3rd pressing of ENTER the following message appears on the display:

CoMeter Calibr Put new sample And press Enter

The buffers' temperature is indicated during measurement. After finishing the measurement the following message appears on the display:

**Co meter Calibr
Put new sample
And press Enter
ADC=xxxxx/1**

Where xxxx is the result from the first calibration measurement. The operator have to put a new buffer, N.B. do not use one and the same buffer more than once! Then start the next measurement. This procedure has to be repeated 5 times. At the end the following message appears on the display:

**CoCalibr-OK
xxxx xxxx xxxx
Xxxx xxxx Diff=xxxx
Power Off/On**

Now the operator has to switch off the power supply of the analyzer. After switching it on again, the analyser has to be cleaned again with water, which ends the calibration of the conductivity measurement system calibration.

Note:

Another possibility for calibration of analyzer's conductivity measurement function.

You need conductivity meter. First measure milk with normal acidity with conductivity meter and use it as sample for calibrating the analysers conductivity measurement function.

4. Co Meter Test.

Serves for testing the working mode of the milk's sample conductivity measurement system. It is used in the production conditions. After this menu is chosen, the analyser executes the procedure for sample's measurement and the display shows the data, used for obtaining the samples conductivity.

**Co Meter Test
CoADC= xxxx
Power Off - Stop**

5. Co Meter En/Dis.

Enables or disables the conductivity measurement system. The following message appears on the display:

Cond Measuring
Yes

No OK
Yes

6. Corrections in conductivity measurement

It is done by starting the menu **Corrections -> Cond measure**. You have the possibility to increase/decrease the measured conductivity value from – 1.00 till +1.00, with step 0.01. After starting this function the display shows the following:

Con Meter
-1.0<=Corr>=1.0
Co Corr=+0.00
Edit – Up/Down

The cursor is positioned below the +. By using buttons **Up/Down**, the operator has the possibility to change the value (number). By pressing the button **Enter**, the operator confirms the chosen value and moves to the next position for editing it. After the last position is edited, if the correction value is within allowed boundaries, the following is displayed: **Co Corr Saved**, which means, that the correction is entered and saved. On the contrary – it returns at the beginning and expects valid correction.

7. Conductivity calibration buffer preparation

In order a standard buffer for conductivity measuring to be prepared follow the instruction below:

1. Take the packet with the powder buffer.
2. Carefully shake the packet in order to gather the powder at the bottom.
3. Cut one end of the packet.

4. Empty its content in a measuring mug with 1 l volume, paying attention all its content to be emptied.

For standard buffer: 5,02 ms – 3,056 r

5. Add 600-700 ml distilled water, which was preliminarily deaerated in vacuum dryer or boiled and then cooled down to 20 °C.
6. Shake the mug till the powder is fully dissolved.
7. Add distilled water to the mark.

Contents:

1. FUNCTION	6
2. TECHNICAL PARAMETERS	7
2.1. WORKING MODES CHARACTERISTICS:	7
2.1.1. MEASUREMENT MODE MILK / DAIRY PRODUCT – FIRST TYPE	7
2.1.2. MEASUREMENT MODE MILK / DAIRY PRODUCT – SECOND TYPE	7
2.1.3. MEASUREMENT MODE MILK / DAIRY PRODUCT – THIRD TYPE	7
2.1.4. CLEANING.....	7
2.2. MEASURING RANGE:.....	8
2.3. ACCURACY:.....	9
2.4 CORRECT AMBIENT CONDITIONS:	9
2.5. DIMENSIONS:	9
2.6. CONTINUOUS WORKING TIME:.....	9
2.7 MILK SAMPLE VOLUME PER ONE MEASUREMENT:.....	9
3. QUALIFICATION OF RAW MILK, THERMALLY TREATED MILK, OTHER DAIRY PRODUCTS AND DERIVATIVES	16
3.1. TAKING SAMPLES AND PREPARATION FOR ANALYSIS	16
3.2. MAKING THE MEASUREMENT.....	16
3.2.1. PREPARING THE ANALYSER FOR WORKING MODE	16
3.2.2. MAKING ANALYSIS	17
3.2.3. DISPLAYING THE RESULTS	18
4. CLEANING THE ANALYSER	21
4.1. AUTOMATIC CLEANING THE ANALYSER	22
4.1.1. CURRENT CLEANING.....	23
4.1.2. END (FINAL) CLEANING.....	25
4.2. PREPARATION OF CLEANING SOLUTION	30
4.2.1. PREPARATION OF ALKALINE CLEANING SOLUTION.....	30
4.2.2. PREPARATION OF ACIDIC CLEANING SOLUTION.....	30
4.3. WORKING WITH THE PERISTALTIC PUMPS	30
5. POSSIBLE MALFUNCTIONS AND ERROR MESSAGES, TROUBLESHOOTING	32
6. MAKING CORRECTIONS AND RECALIBRATION OF THE DEVICE.....	33
6.1. TAKING SAMPLES AND PREPARATION OF SAMPLES FOR CHECKING THE ACCURACY OF THE MILK ANALYSER, MAKING CORRECTIONS AND RECALIBRATION	33
6.2. DETERMINATION THE TYPE OF THE DISCREPANCY:	33
6.2.1. MAKING MEASUREMENTS	33
6.2.2. ANALYSING THE MEASUREMENT RESULTS.....	33
6.3. MAKING CORRECTIONS	34
6.3.1. POSSIBLE CORRECTIONS, LIMITS AND CHANGING STEPS.....	34

6.3.2. PREPARING THE ANALYZER FOR MODE CORRECTIONS	35
6.3.3. MAKING CORRECTION.....	35
CHECKING THE CALIBRATION	42
7. STARTING THE DEVICE IN A SERVICE TEST/SETUP OPERATIONAL MODE. MENUS DESIGNATION.....	43
7.1. STARTING THE DEVICE IN A SERVICE TEST/SETUP OPERATIONAL MODE.....	43
7.2. MENUS FUNCTION:	44
7.2.1. SPECIAL MODES.	44
7.2.2. CORRECTIONS	44
7.2.3. SETTINGS.	44
7.2.4. TESTS.....	50
7.2.5. PH METER & CO METER	51
7.2.6. EXIT	51
7.2.7. SUMMARY OF MILK ANALYSER SETUP FUNCTIONS	52
8. ADDITIONAL POSSIBILITIES OF THE ANALYSER.....	55
8.1. CONNECTING TO 12 V DC POWER SUPPLY.....	55
8.2. CONNECTING TO IBM PC.....	55
8.3. CONNECTING ADDITIONAL SERIAL PRINTER (OPTION).....	55
8.4. CONNECTING MATRIX PRINTER (OPTION).....	55
8.5. CONNECTING AND WORKING WITH EXTERNAL KEYPAD (OPTION).	55
8.6. MEASURING HIGH FAT SAMPLES (OPTION).....	57
8.7. WORKING WITH EMBEDDED ACCUMULATOR (OPTION).....	58
8.8. EMBEDDED REAL TIME CLOCK (OPTION).	59
9.1.1. USB FLASH-DRIVE (PEN DRIVE) (OPTION).	60
APPENDICES	62
APPENDIX 1: PREPARATION OF SAMPLES FOR MILKANALYSERS' CALIBRATION.....	62
1. DETERMINATION SOME OF THE PARAMETERS BY FORMULAS.....	64
2. SNF DETERMINATION.	64
3. DETERMINATION OF LACTOSE CONTENT	64
4. DETERMINATION OF SALTS CONTENT.....	65
5. DETERMINATION OF TOTAL PROTEINS CONTENT	65
APPENDIX 2 FREEZING POINT DETERMINATION	67
1. METHODS FOR DETERMINATION.	67
2. THE BASIC FREEZING POINT.	67
APPENDIX 3 PH MEASURING (ONLY FOR DEVICES WITH EMBEDDED FUNCTION).....	70
1. GENERAL INFORMATION	70
2. PH ELECTRODE.....	70
2.1. ELECTRODE PART	70
2.2. ELECTRODE CARE & ELECTRODE MAINTENANCE	70
2.3. STORAGE.....	70
2.4. AFTER USE.....	70

2.5. ELECTROLYTE REPLACEMENT (FOR REFILLABLE ELECTRODE ONLY).....	71
2.6. NEW ELECTROLYTE PREPARATION:	71
2.7. RE-USE THE ELECTRODE.	71
2.8. ELECTRODE CLEANING.....	71
2.9. ELECTRODE ACTIVATION	72
2.10. REJUVENATION PROCEDURE	72
2.11. ELECTRODE LIFESPAN	72
3. BUFFER SOLUTIONS	73
4. PH ELECTRODE CALIBRATION	73
5. PH HELPFUL HINTS	74
1. PREPARATION FOR PH MEASUREMENT	76
FIG. 11 PLACING THE STOPPER.....	76
FIG. 12 PLACING THE PROBE.....	77
7. PH MEASURING.....	78
8. <i>pH Calibration.</i>	80
FIG. 13 CALIBRATION BUFFERS	80
9. <i>pH Measuring.</i>	81
10. <i>pH En/Disable</i>	82
11. <i>pH U Display</i>	82
12. <i>pH test</i>	82
APPENDIX 4 CONDUCTIVITY MEASURING (ONLY FOR DEVICES WITH EMBEDDED FUNCTION)...	83
1. METHOD OF DETERMINATION.....	83
2. CONDUCTIVITY MEASUREMENT.....	84
3. CO METER CALIBR	84
4. CO METER TEST.	85
5. CO METER EN/DIS.	86
6.CORRECTIONS IN CONDUCTIVITY MEASUREMENT	86
7. CONDUCTIVITY CALIBRATION BUFFER PREPARATION	86
CONTENTS:.....	88

GUARANTEE CARD

LACTOSCAN LA

**Guarantee period is 1 (one) year after purchasing date.
Improper handling, transport and storage will invalidate the guarantee.
Guarantee is void if warranty labels are removed.**

Serial №

Date of purchase:

Password:

Distributor:

Signature:

Stamp:

GUARANTEE CARD**Purchaser:**

--

Service report:

Service entry date	Damage	Delivery date	Signature

Covers:**Lactoscan LA****Last edited: 30.08.2016****HEAD OFFICE**

4, Narodni Buditeli Str.

8900 Nova Zagora, BULGARIA

Phone/Fax: + 359 457 67082

e-mail: office@lactoscan.com

www.lactoscan.comwww.milkotronic.com<https://www.facebook.com/lactoscan><https://www.youtube.com/user/LACTOSCAN><https://plus.google.com/+LACTOSCAN><https://www.linkedin.com/company/lactoscan>