

Lactoscan TETA

Thermostatic device with **AUTO-DROP** System

Simple and easy to use device with auto-drop system for incubation of tests for detection of antibiotics in milk.

- 2 working modes:**
- **automatic (for 5 tests)**
 - **manual (for 14 tests)**

KEY FEATURES:

- **User-friendly: simple in operation, setting, maintenance and installation.**
- **Portable and compact design**
- **Stainless-steel, shockproof body**
- **Low power consumption**
12V/25W and 24V/100W (only during the heat-up process)
- **Mobile (12 to 24 V DC power supply)**
- **Easy cleaning (stainless-steel body)**
- **Light weight (only 0.8 kg)**
- **OEM (laser engraved branding)**

The strip holder has two positions:



1st when placing the test vials into the thermostatic device.



2nd when placing the test strips into the strip holder.

Specification:

Model name	Number of wells	Compatible test/tube:
Lactoscan TETA	5/14	Unisensor's milk tests: Twinsensor, 4Sensor, AflaSensor Quanti, BetaXpress, CapSensor, GoatSensor, MelamineSensor, QuinoSensor, SulfaSenSor, TyloSensor; Ballya Bio's milk tests; BioEasy's milk tests; Zeulab's milk tests: Eclipse 3G .

Temperature range:
30°C to 80°C (± 0,5°C)

Time range:
• 0 to 244 minutes



**ENVIRONMENTAL
CONDITIONS:**

Ambient air temperature
has to be 5°C below
the desired set temperature

Operating temperature:
-15°C to 45°C

ELECTRICAL PARAMETERS:

Switching adapter:

Input:
100 - 240 V ~ 1.6 A max
50 - 60 Hz

Power consumption:
12V/25W and 24V/100W
(only during the
heat-up process)

**MECHANICAL
PARAMETERS:**

DIMENSIONS
(W x L x H):
manual mode: 115 x 140 x 45
automatic: 115 x 150 x 90

Weight: 0,8 kg

Stainless steel
cover box

buy online at
www.lactoscan.com



Lactoscan TETA


Thermostatic device with **AUTO-DROP** System

High accuracy and easy to use

Portable and programmable

Uniform and stable temperature

Thermostatic device can be
connected towards **truck battery 12V - 24V**

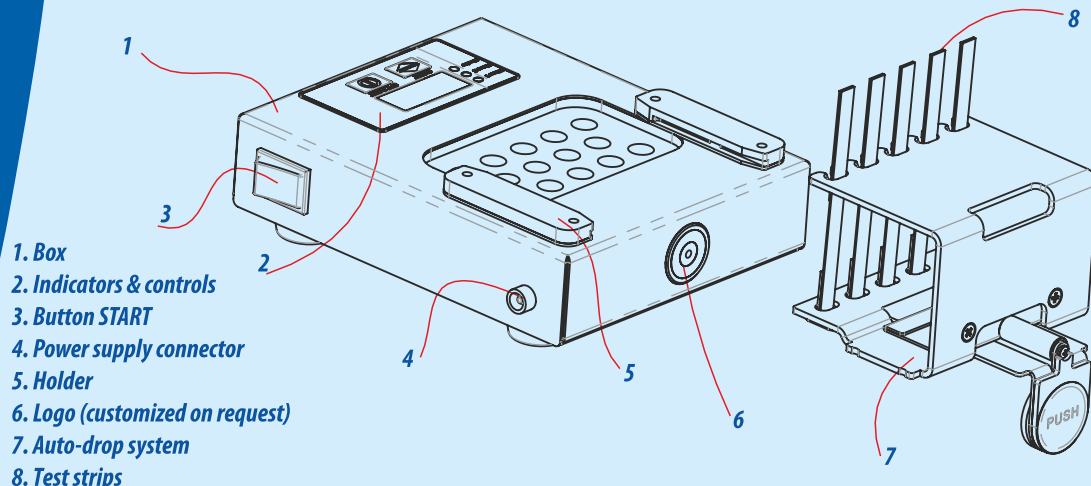

On customers' request **THERMOSTATIC DEVICE**
can be **OEM (laser engraved branding)**



MILKOTRONIC LTD

Bulgaria, 6000 Stara Zagora
167, Tsar Simeon Veliki, A, 25
Tel: + 359 457 670 82,
Fax: + 359 457 670 83

e-mail: salesdpt@lactoscan.com
Eng, Spa, Port: +359 886 841 389
e-mail: sales2@lactoscan.com
Eng, Rus: +359 886 841 395
e-mail: mainservice@lactoscan.com
Eng, Rus: +359 888 811 017

e-mail: m.mihneva@lactoscan.com
Eng: + 359 885 597 232
e-mail: milkanalyzers@lactoscan.com
Eng, Rus: +359 896 753 803

distributors: